

# ALGINATES, SPATULAS, BOWLS & KNIVES

## Alginmax

THE ONLY CHROMATIC-PHASE HYDROCOLLOID ALGINATE FOR ORTHODONTICS, CROWN / BRIDGE AND DENTURES

### HIGH PRECISION ALGINATE FOR IMPRESSIONS WITH CHROMATIC PHASE INDICATOR

Thirty years ago we were the first to introduce a chromatic alginate to the United States and now in our European facility, we are manufacturing an alginate designed for dentures, crowns, bridges and orthodontics, with these unique features:

- **Dimensional Stability** - more than 120 hours (After washing, rinsing and removing excess water by shaking, impression must be placed in a zip-loc bag.)
  - **Pleasant Scent** - Patients like the pleasant vanilla scent which has a neutral taste.
  - **Chromatic Phases:**
    - **Violet** (mixing - approx. 35 seconds)
    - **Rose** (packing/working - approx. 40 sec.)
    - **Light Blue** (setting - approx. 45 seconds)
- The impression is complete in just about 2 minutes**
- Final impression is light blue for easier reading
  - Lead and cadmium free
  - Very low dust
  - Bio-compatible
  - Meets ISO 1563 and ADA #18
  - Easy removal - high tear strength
  - Four pound package includes:
    - 4-one lb. bags and 2 measuring cups
  - 510k certificate on file



GROUP/ITEM #	DESCRIPTION
003-2802/1	Alginmax Container w/1 lb. Bag + Cups
003-2802/4	Alginmax Box of 4 -1 lb. Bags
003-2802/S2	Alginmax 100 gr. sample w/ set of cups
003-2802/CUPS	Vials for Powder & Liquid (1 set)
003-2803	Alginmax Storage Container



## Nylon Spatulas for Alginates

GROUP/ITEM #	DESCRIPTION
003-2803/BS	Alginmax Blue Spatula
003-2803/WS	Alginmax White Spatula



## Spatulas for Alginate and Plaster

This line of spatulas was designed by dental engineers for better manipulation of alginate and plaster. The curved tips on items 017-434, 017-435, and 017-436 are the most advanced technical design. They help to mix and become a scoop which is a great help in picking up additional material.

GROUP/ITEM #	DESCRIPTION
017-430	7 1/4" Spatula
017-431	7" Spatula
017-432	8 1/4" Spatula
017-433	8" Flexible Spatula
017-434	7 3/4" Flexible Spatula
017-435	7 3/4" Rigid, Heavy Blade Spatula
017-436	8 3/4" Rigid, Heavy Blade Spatula



017-430

017-431

# Pozzi®



Since 1971

017-432



017-433

017-434

017-435

017-436

# Pozzi Super-Flexible Mixing Bowls

Their soft flexible texture and sturdy construction make these bowls indispensable for mixing plaster or alginate. They easily conform to the hand and make mixing a smooth and quick process. These self cleaning plasticized polyvinyl bowls are available in four convenient sizes from 350ml to the max 2250ml size.

GROUP/ITEM #	DESCRIPTION
017-5501-10	Flexible mixing bowl small blue
017-5501-12	Flexible mixing bowl regular blue
017-5501-14	Flexible mixing bowl large blue
017-5501-16	Flexible mixing bowl max blue




**Pozzi**<sup>®</sup>  
Since 1971

## Plaster Knives

High quality and fine craftsmanship set these knives apart.

GROUP/ITEM #	DESCRIPTION
017-426	Plaster Knife, 6 3/4"
017-428	Plaster Knife, 6 1/4"
017-429	Plaster Knife, 6 1/4"

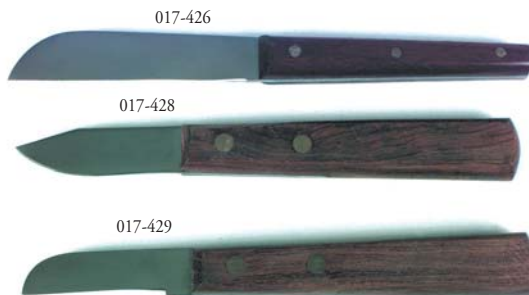

## Special Inox Plaster Nippers

Stainless steel nippers especially handy for accurate chipping of excess plaster from laboratory work models. Cuts plaster and stone.

GROUP/ITEM #	DESCRIPTION
007-444	Special Inox Plaster Nippers



## Laboratory Knives



GROUP/ITEM #	DESCRIPTION
017-420	Fahnstock Large Waxing Knife 7"
017-421	Fahnstock Small Waxing Knife 5 1/4"
017-422	Lessmann Small Waxing Knife 5 1/4"
017-423	Lessmann Large Waxing Knife 7"
017-424	Waxing Knife with Large Spoon 7"
017-425	Waxing Knife with Small Spoon 5"

## Flask Opening Knife

The dual purpose flask opening knife is a heavy duty, long lasting instrument made of carbon steel that extends the length of the handle. The butt end has an extrusion designed to open flasks.

GROUP/ITEM #	DESCRIPTION
017-427	Flask Opening Knife 6 3/4"



## Knife & Spatula Kit

A complete kit of 18 different high quality instruments, wrapped in a protective red felt roll.

GROUP/ITEM #	DESCRIPTION
017-4500	Knife and Spatula Kit



**Pozzi**  
Since 1971