

A Brand Division of American Tooth Industries

1200 Stellar Drive • Oxnard, CA 93033-2404
(805) 487-9868 • Toll Free (800) 235-4639 • FAX: (805) 483-8482
E-mail: info@americantooth.com • Web Site: www.americantooth.com



AVM 100M ARTICULATOR


Characteristics and Uses.

The AVM 100M Articulator from Major America is a medium-priced articulator with no face-bow and it is extremely easy to use. Simulation of lateral and protrusive excursions is simplified since the models are mounted using arbitrary standard values. It can be used for all traditional applications including full dentures, crown and bridge work and partials.

The AVM 100M Articulator is extremely light-weight making it very manageable while in use. Maximum aperture angle is 120 degrees. The AVM 100M allows for complete visibility while working on pieces.

Dimensions and Technical Characteristics:

- Height: 95 mm
- Intercondylar width: 105 mm
- Bennet angle set at 23 degrees
- Includes mounting plane and three mounting plates
- Incisal table reference punch market at zero
- Parallelism of the instrument is obtained by setting the instrument at the punch mark.

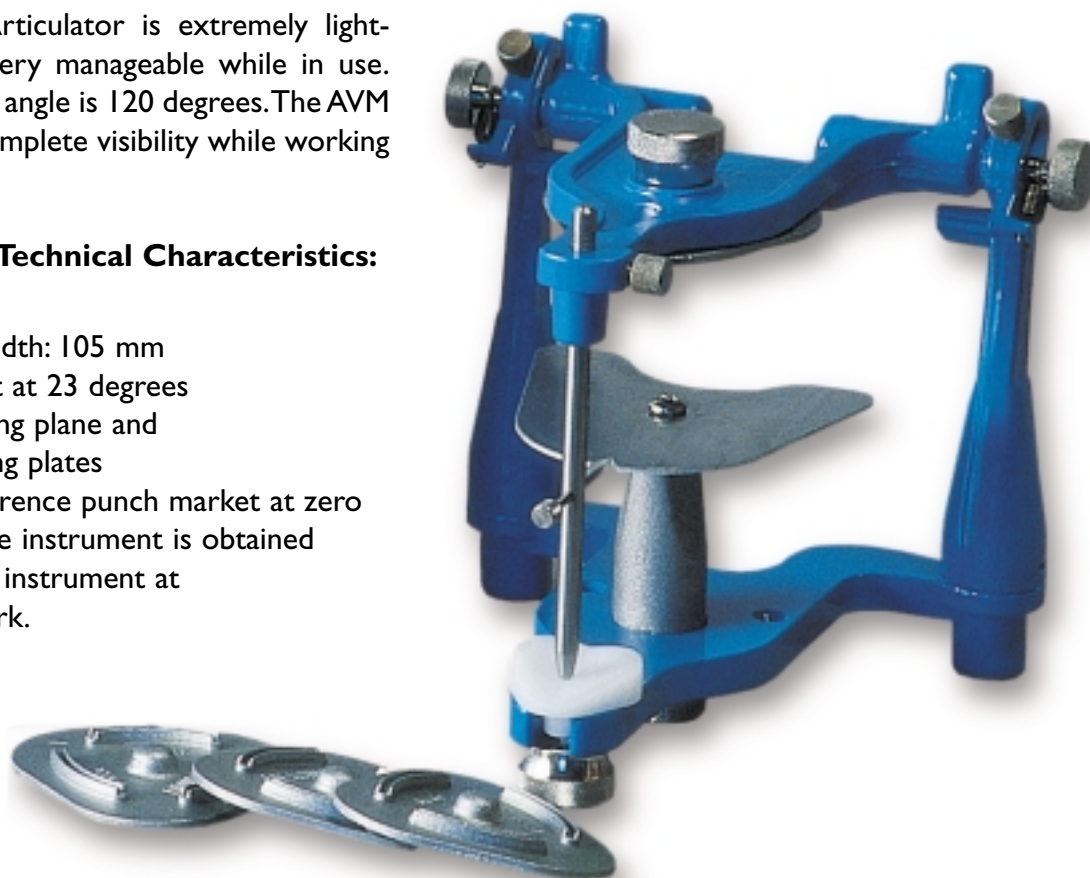

PACKAGING AVAILABLE:

Grp/Item #	Description	Quantity
073-E3000	AVM 100M	One Unit

**Additional questions, to place your order, or for the Dealer nearest you, please call the
Tooth Express Order Desk at 800-628-1437.**

