

# AMERICAN MADE TEETH LARGEST TOOTH INVENTORY AVAILABLE IN NORTH AMERICA



### **Table of Contents**

CAD/CAM Material	3-5
Justi Teeth Products	6-17
Justi Specialty Teeth Products	18-20
Justi Tooth Cabinets and Tray Specials	21-22
Justi Denture Resin	23
Acrylic Materials for Prosthetic use	24-25
Acrylic Materials for Dental Lab use	26-27
Resin for Prosthetic use	28
Specialty Resin and Accessories	29
Resin for Orthodontic use	30-31
Accessories for Denture Maintenance	32
Accessories for Making Dentures	33
Alginates for Impressions	34-35
Accessories for Making Impressions w/Silicones	35
Condensation Curing Silicones	36-39
Wax	39
Pozzi Super Tray Systems/SS Impression Trays	40
Pozzi Impression Trays and Accessories	41
Laboratory Specials	42-43
Dental Porclein Products	44
Galletti Articulator	45
Articulators	46
Equipment and Laboratory Specials	47
Microwaveable Flasks	48
Pozzi Highest Quality Trimmers	49
Pozzi Abrasive Universal Series	50-56
Pozzi Stainless Steel Crowns	57-58
Pozzi Straight Length and Spool Wires	59
Pozzi Wire and Silver Solder	60-61
Pozzi Utility Pliers	62-70
TCS Dental Supply Products	71-73
ASA Dental Products	74
Galloni Crucibles	75
DentalFarm	76-79



**Justi's Top Line!** For quality, durability & price, they're the best teeth in the world! Many universities use them as a test standard. Recommended for your "Premium" House Brand Cases. Premium quality at a meduim price!

Umperial VIVO	AS LOW AS	
	\$16.02	\$14.31
	ANTERIOR 1x6	POSTERIOR 1x8

Justi Imperial ViVO<sup>®</sup> offers traditional quality American molds using the newest Vita shade system. The effect of this combination is outstanding aesthetics and excellent functionality.

DYMON-HUE HPT 2N Dymon-Hue HPT H58		AS LOW AS	
	291 Dymon-Hue HPT HIS2	\$6.05 ANTERIOR 1x6	<b>\$6.05</b> POSTERIOR 1x8
	A25 Dymon-Hue V A2	AS LO	W AS
DYMON-HUE	32M Dymon-Hue V A2	\$3.08	\$3.08
•	0.0444400	ANTERIOR 1x6	POSTERIOR 1x8

**Quality at an Economy Price!** For the price concerned labs that want to maximize quality for their "House Brand" teeth. Available in Bioform and Vita shades. Quality at an economy price!

JUSTI BLEND®	AS LO	AS LOW AS	
	\$18.92	\$14.31	
0011100	ANTERIOR 1x6	POSTERIOR 1x8	

For the lover of the natural look. Justi Blend molds are similar to Bioblend, Trublend, and Portrait. Recommended for all your "Blended Premium" cases. Premium quality at a medium price!

### **CAD/CAM Material - A Complete System**

### Justi® High-Impact Denture Base PMMA Discs





### Justi High-Impact Denture Base PMMA Discs

are intended for use in milling full and partial denture prosthetic appliances. Justi Discs are available in 98.5mm diameter to fit most open type, 5-axis milling machines. Discs are available in three heights, 20mm, 25mm, and 30mm and are available in four shades, light pink, light red pink, dark pink, and regular original.



# "Outstanding aesthetics and excellent functionality."

- Available in 19 of the newest Vita<sup>®</sup> type shades, including 3 bleached
- Available in 24 upper and 9 lower anterior molds, 17 upper and 17 lower posterior molds,
- Made in the U.S.A. backed by over 150 years of experience



See page 10-11 for specials

### Justi<sup>®</sup> Bonding/Repair Material

- Universal type bonding/repair material
- Self-curing bonding/repair methyl methacaylate
- · Cross-link and internally pigmented to give a dense even texture
- Available in 3 shades of powder: regular pink, regular fiber pink and clear

We recommend using the Justi Bonding/Repair Material as it is compatible with the Justi Denture Base PMMA Discs and Justi Imperial ViVO Denture Teeth



#### See page 24 for pricing & specials

# A complete system that can't be beat

Jusit Imperial ViVO Digital Denture System: by using Justi's denture material together, Justi produces the best posible final appliance. Justi's Dentue Base PMMA Disc and Justi Imperial ViVO are bonded together using Justi Repair Material provides a dense even texture to the denture teeth.



call today to place your order 805-487-9868

# Flexible Clear Disc 🗲

For Invisible Partial Denture

The Flexible Disc is intended for aesthetically superior semi-flexible removable partial denture frameworks, which blend perfectly with the patient's gums, providing better function, aesthetics, comfort, and durability.



### **Benefits for Labs:**

- 100% digital for accurate fit
- Eliminates RPD in the wax-up stage, investing, in casting and finishing.
- Fast, Clean, precise, and cost-effective
- Easy finishing and polishing with minimal effort
- Mill with the same procedure with PMMA, both wet and dry mill

### Indications:

Transparent partial denture frameworks, night guards, clasps, and other dental appliances processed by CAD/CAM milling machines for lightweight, aesthetic nonmetallic restorations.



### **Benefits to Patients:**

- Lightweight, durable and precise
- Comfortable and easy adaptation in the mouth
- Bio-campatible and non-allergenic
- Clasps are flexible and tooth-bearing, virtually unbreakable.
- Occlusal rest reduce tissue-borne pressure to reduce sore spot and offer excellent retention

GROUP/ITEM #	DESCRIPTION	
102-FLD9820-CL	20mm x 98.5 mm (2 steps) – Shade: Transparent	\$70.00

### **Full Denture PMMA Disc**

More than a temporary material.

- Reduce production process.
- Denture Arch is stable on the denture base.
- Better accuracy than hand-made full denture.

### Shades Available:


Denture Base: 2S, 3ST, 3S, 3ST Denture Arch: A1, A2, A3, A3.5, B1, B2, B3, C2, D2

GROUP/ITEM #	DESCRIPTION	PRICE
102-2S9820/Shade*	Full Denture PMMA Disc 98X20	\$92.00
102-2ST9820/Shade	Full Denture PMMA Disc 98x20	\$92.00
102-3S9820/Shade	Full Denture PMMA Disc 98x20	\$92.00
102-3ST9820/Shade	Full Denture PMMA Disc 98x20	\$92.00
102-2S9825/Shade	Full Denture PMMA Disc 98X25	\$95.83
102-2ST9825/Shade	Full Denture PMMA Disc 98x25	\$95.83
102-3S9825/Shade	Full Denture PMMA Disc 98x25	\$95.83
102-3ST925/Shade	Full Denture PMMA Disc 98x25	\$95.83
102-2S9830/Shade	Full Denture PMMA Disc 98X30	\$101.58
102-2ST9830/Shade	Full Denture PMMA Disc 98x30	\$101.58
102-3S9830/Shade	Full Denture PMMA Disc 98x30	\$101.58
102-3ST9830/Shade	Full Denture PMMA Disc 98x30	\$101.58

\*example 102-2S9820/A1





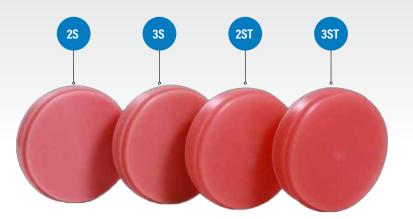

### **CAD/CAM** Material

### Justi<sup>®</sup> High-Impact Denture Base PMMA Discs





\$52.88 30mm and are available in four shades, light pink, light red DB PMMA Disc 98mm x 30mm \$58.18 DB PMMA Shade Guide \$8.50 pink, dark pink, and regular original.

\*Specify shade when ordering –example 103 DB20-DKPK

### Justi<sup>®</sup> Monolayer PMMA Disc

Made from high cross linked PMMA material, the monolayer PMMA disc are used for fully and partially anatomical long-temporary crowns and bridges.

#### Advantages:

103-DB30-Shade\*

103-DBPMMA-SG

- Biocompatibility
- Excellent polish ability
- Clean milling
- Excellent long-term shade stability and aesthetic
- Available in shades: Pink 2S, A0 White, A1, A2, A3, A3.5, A4, B1, B2, B3, B4, C1, C2, C3, C4, D2, D3, D4



GROUP/ITEM #	DESCRIPTION	
102-HADB35-shade	Justi PMMA Disc 98x14mm	\$43.70
102-HADB39-shade	Justi PMMA Disc 98x20mm	\$43.70
102-HADB41-shade	Justi PMMA Disc 98x25mm	\$49.45
*Specify shade when o	ordering –example 102-HADB35-2S	

# **10 LAYERS**

### Justi<sup>®</sup> Multilayer PMMA Disc

Fully adjustable partially anatomical long-term temporary crowns and bridges, especially for the anterior area.

#### **ADVANTAGES**

- Natural tooth-like layering 10 layers
- For long-term temporaries with a wear period up to 12 months
- Biocompatibility
- Shade stability and plaque build-up
- Excellent polishing performance after milling

SHADES: A1, A2, A3, A3.5, B1, B2, B3, C2, D2, BL (Bleached Shades)



		Special Intro Price
GROUP/ITEM #	DESCRIPTION	
102-HAMDB21-Shade	Justi Multilayer PMMA 98x14mm	\$80.50
102-HAMDB25-Shade	Justi Multilayer PMMA 98x20mm	\$80.50
102-HAMDB27-Shade	Justi Multilayer PMMA 98x25mm	\$86.25
*Specify shade when ordering example HAMDB21-A1		

For Beautiful Smiles

# **DYMON-HUE® HPT**

<section-header>

The posterior line is available in multiple molds within the 0°, 10°, 20° and 33° occlusal forms with functional Justi<sup>®</sup> design.

# **DYMON-HUE**<sup>®</sup> HPT

Anteriors and Posteriors Super Quality at an Economy Price



### **DYMON-HUE**<sup>®</sup> **HPT** Single Teeth

For added savings, laboratories can purchase uncarded single teeth. Minimum quantity is 500 teeth per order. Please allow ten days delivery time for orders of single teeth.

 GROUP/ITEM #
 DESCRIPTION

 081-19973-1
 Dymon-Hue HPT Anterior \$1.67

 081-19974-1
 Dymon-Hue HPT Posterior \$1.24



6



### Justi<sup>®</sup> Dymon-Hue<sup>®</sup> Hardened Plastic Tooth line now includes Extra Large Anterior Mold #268.

**DYMON-HUE HPT** teeth are manufactured with life-like aesthetics that are blended to be similar to those of natural teeth.

**Dymon-Hue HPT** (hardened plastic teeth) are made of high quality synthetic resin. They are wear resistant, cross-linked, vacuum processed, and have superior hardness. The cross-linking system used is an interpenetrated, multi-layered, intersectional polymer network that is solvent and heat resistant. The pigmented copolymer incorporates inorganic type pigments that are color stable and unaffected by light or heat. The denture base side of the tooth extending to the neck portion is constructed of a solvent- resistant, softer polymer that assures a 100% chemical bond to all methyl methacrylate denture base acrylics. No diatoric retention is necessary.

**Dymon-Hue HPT** teeth are available in 35 anterior molds, and 40 posterior molds including 0° Planatomical<sup>™</sup>, 10° Symmetry<sup>™</sup>, 20° Semi Anatomical, and 33° Anatomical and Short-Bite. They are also available in anatomical short bite molds that provide the technician additional space, eliminating tooth adjustments in many cases.

The **Dymon-Hue HPT** is another quality Justi<sup>®</sup> product manufactured to surpass ANSI/ADA Specification #15 for synthetic resin teeth; it also complies with ISO # 22112:2005.



Extra Large Anterior Mold Available: #268

# Available as a FREE download at www.americantooth.com/dhmc



HOUSE BRAND PROMOTION-NEW AND DORMANT ACCOUNTS

## **Justi® Teeth Products**

A2

For Beautiful Smiles





# **DYMON-HUE<sup>®</sup> HPT New Bleached Shades Dymon-Hue V** Anterior/Posterior **NOW AVAILABLE** Dymon-Hue<sup>®</sup> Maxillary Centrals Shade Guide Front

Two new striking bleached shades have been added to the Dymon-Hue Tooth line & Shade Guide: Shades H21 & H22.

Group/Item # Description Dymon-Hue HPT SG 15 Shades Double-sided \$24.40 084-19987



Dymon-Hue is available in Vita\* Shades & 12 popular molds.

32M

Group/Item # 085-19973/V 085-19974/V

Description Dymon-Hue V Anterior 1x6 \$3.85 Dymon-Hue V Posterior 1x8 \$3.85

### Dymon-Hue V 12 Tab Shade Guide

Dymon-Hue V tooth line is available in 12 Vita\* Shades:

Back

A0, A1, A2, A3, A3.5, B1, B2, B3, C1, C2, C3, D2

Group/Item # Description 085-19975/SG Dymon-Hue V 12 Tab Shade Guide \$3.50 \*Vita Shades' is not a trademark of American Tooth Industries.

A35 **B1 B2 C**3 D2 A0 A2 A3 **B**3 C1 C2 A1 Dymon-Hue 

Call 805-487-9868 to request a free shade guide

For Beautiful Smiles



Justi<sup>®</sup>, the professional's choice for hardened plastic teeth, celebrates over 158 years in the dental industry!



- Clean Flash Lines Less Grinding
- Vacuum Processed Denser Acrylic
- Excellent Bond Strength Fewer Repairs
- Exceptional Aesthetics Patient Satisfaction

GROUP/ITEM #	DESCRIPTION	
084-51-0425	Imperial Age 18 Shades	\$29.28
084-51-0426	Imperial Age Combo Shade Guide Set 18+4 Shades	\$35.80
084-51-0425TS	Imperial Bleached 4 Shades	\$6.52

HOUSE BRAND PROMOTION-NEW ACCOUNTS

**22 SHADES** Based on Age **INCLUDING NEW BLEACHED SHADES** 

USD American Mode Teeth COLOR QUALITY CONTROLLED BY

🗙 x-rite PANTONE®

lusti

57 Anterior Molds - 52 Posterior Molds



Imperial® Age Shades



Imperial<sup>®</sup> Bleached Shades



Since 1864

## **IMPERIAL NEW USER PACKAGE**

# **IMPERIAL MAIN FEATURES**



**SEE PAGE 20 FOR DISCOUNTS** 


805-487-9868




For Beautiful Smiles

	Ömperial Special	Promotion	S	
-	20 H.D. JUSTI CO.	Cashing		
BUY	10 Cards Imperial Teeth	BUY	\$1,207.50 of Imperial Teeth	
GET	1 CARD FREE	GET	<b>GET</b> FREE JUSTI TIMELESS ARTICULATOR	

Delta  $\Delta$  - The Next Highest Quality Generation of Imperial Molds

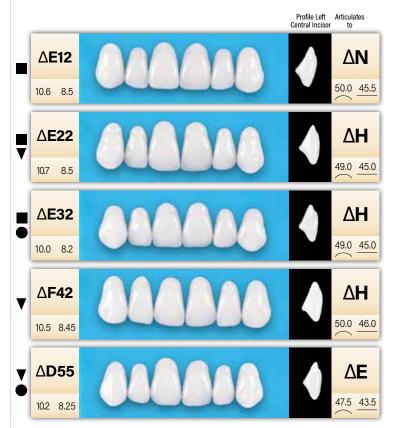

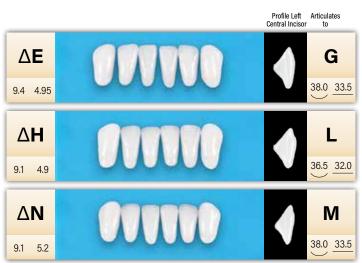

### For articulating posteriors see page 16

For Beautiful Smiles



*Justi*°, the professional's choice for hardened plastic teeth, celebrates over 158 years in the dental industry!



"Outstanding aesthetics and excellent functionality."



Justi Imperial ViVO<sup>®</sup> offers an incisal and interproximal blend that give this beautiful line of teeth the unique ability to absorb and reflect surrounding color. This means more latitude for blending with natural teeth resulting in less resets due to shade, especially in partial denture cases.

Justi Imperial ViVO<sup>®</sup> offers traditional quality American molds using the newest Vita shade system. The effect of this combination is outstanding aesthetics and excellent functionality.


The ability of the teeth to blend into the adjacent dentition has made producing great-looking stay plates much more predictable.

The chameleon-like characteristics of Justi Imperial ViVO<sup>®</sup> teeth make shade selection much easier and greatly reduces remakes. The traditional American molds have a very natural look to them.

### Imperial ViVO<sup>®</sup> teeth are suited for:

- Full dentures
- Implant prosthetics
- Telescope and cone constructions
- Attachment dentures
- Model cast prosthesis
- Interim solutions

# Delta $\Delta$ - The Next Highest Quality Generation of Imperial ViVO Molds

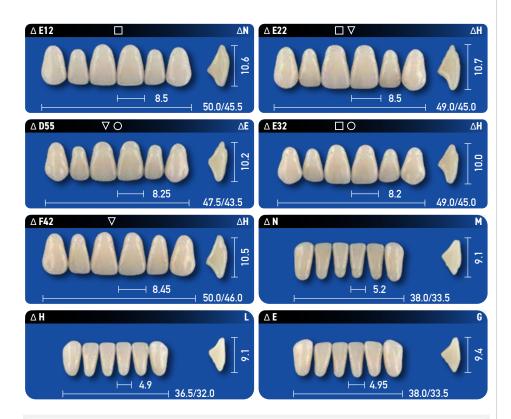

### For articulating posteriors see page 16

### IMPERIAL ViVO<sup>®</sup> Specialty Tooth Set Items

### **Upper Anterior and Posterior Extra Large Teeth**

Imperial ViVO<sup>®</sup> Special upper anterior AXL central incisors and upper posterior PXL36 bicuspids are specialty items that will save you the expense and trouble of making custom teeth when you need to fill a large anterior space in a complete upper or partial denture. The upper posteriors may also be used as molars. These specialty teeth are available in all 19 ViVO<sup>®</sup> shades.

BUY	3 Tooth Cards Central Incisors
GET	1 BICUSPID FREE

### **Ideal Flipper Saver**



Made in the USA | All teeth are made in the USA backed by over 158 years of experience. All raw materials purchased from US factories.

10

FOR REORDERS SEE PAGE 20 FOR DISCOUNTS

For Beautiful Smiles



lusti

COLOR QUALITY CONTROLLED BY

xrite PANTONE®

805-487-9868

For Beautiful Smiles

**BUY** 

**GET** 

**GET** 









### Justi<sup>®</sup> Blend Mold Chart



Tooth Cards (1x6 or 1x8)

HOUSE BRAND PROMOTION-NEW AND DORMANT ACCOUNTS

Justi<sup>®</sup> Blend New User Package

**1 FREE TOOTH CABINET** (with Inserts)

# Justi<sup>®</sup> Blend Stocking Special



26 Posterior Cards FREE

# **Justi<sup>®</sup> Blend MAIN FEATURES**

The structure of Justi Blend<sup>®</sup> anterior hardened plastic teeth creates the vitality of natural teeth. The centrals, laterals and cuspids are manufactured with the same HUE, but each type has a different VALUE. It is this combination of tones blended together in one set that gives that natural look and natural characterizations.

The semi-translucent body creates the look of depth within the tooth and reflects color from within, as dentine does in a natural tooth.



The Justi Blend<sup>®</sup> Posterior visually matches as closely as possible the cuspids, left and right. It does not match the value of the centrals or the lateral teeth. It is manufactured with the same HUE but has different VALUES or TONES. Following the studies of U.S. universities, it was concluded that natural posteriors are sightly darker than the centrals and laterals. American Tooth Industries is following this characteristic.

Offer Expires December 31, 2024–not valid with other offers www.americantooth.com | buy online: store.americantooth.com

Justi<sup>®</sup> Teeth Products For Beautiful Smiles

Made in Europe

**SPECIAL** 

Buy Super-Lux & Major Dent Anteriors at \$18.76 and get Triplex Posteriors at \$13.50

# SPECIAL FEMALE & MALE MOLDS ABSOLUTELY UNIQUE! Triplex™ MSLV

### The European Look in Vita type of shades

The Triplex<sup>™</sup> MSLV is a unique system of multilayered intersectional polymer network resin denture teeth. The Triplex System features the natural-looking characterization found in European manufactured female and male anteriors, along with the full range of anatomical occlusion in the posteriors required by modern dental professionals in the USA. The Triplex System is a unique line of teeth for modern users in over-denture cases or partial cases where high characterization is needed. It is also the only tooth line in the world market to offer male and female anterior maxillary and mandibular molds for very exclusive premium full dentures.


#### Available in 12 shades A1, A2, A3, A3.5, A4, B2, B3, B4, C1, C2, C4, & D4

GROUP/ITEM #	DESCRIPTION	
085-01477/V	Super Lux V Anterior Card (1x6)	\$20.84
085-01577	Major Dent V Anterior Card (1x6)	\$20.84
082-51-0101	Justi Triplex MSLV Posterior Card (1x8)	\$17.89
084-01578/SMT	Shade Guide (for use with the Triplex System) 12 Vita Shades	\$19.53

# TRIPLEX MSLV MAIN FEATURES

A unique system featuring European aesthetic and American technology.

Major Dent<sup>®</sup> female anterior molds designed and sculptured in Italy.

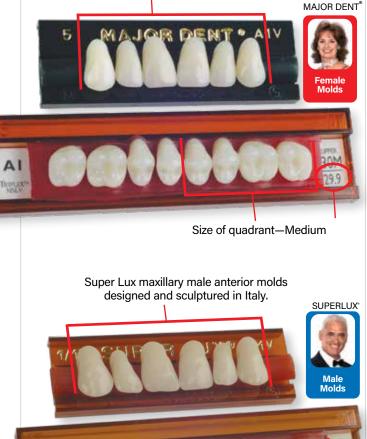



Size of quadrant—Medium

Triplex posteriors are designed and manufactured in the United States and are available in 0°, 10°, 20°, and 33° to articulate to the Super Lux and Major Dent anteriors. The Triplex<sup>™</sup> MSLV System utilizes a Vitatype shade system.

Offer Expires December 31, 2024-not valid with other offers www.americantooth.com | buy online: store.americantooth.com



Made in Four Lavers, Dentine

Acrylic and Enamel

### SPECIALIZED COMPOSITE TEETH

# Major<sup>®</sup> Plus Comp

### Anterior Molds and Posterior Molds

Major Plus Comp teeth with enamel layer in composite material; dentine and cervical layers in PMMA acrylic materials to ensure a perfect chemical bonding to the denture base.

- Very high surface hardness
- Excellent resistance to abrasion
- Precise occlusal profiles
- Assortment of natural molds and shades
- Permanent color stability
- Cuspidal inclination of Anatomical molds between 23° and 25°

RETAIL

PRICE

\$ 32.21

\$ 32.21

\$ 32.21

\$ 55.84

SPECIAL

PRICE

\$21.45

\$21.45

\$21.45

\$48.91

- Cuspidal inclination of Technical molds are around 20°
- 4 layers acrylic and composite teeth



Easier Bondability With Tubular Neck Design

וס P MAIOR 133 Shade Number Posterior Molds come in three occlusal profiles: High Anatomical, Technical, and Wide Functional Area Major Plus Comp Anterior teeth are offered in 10 upper and 4 lower molds in 3 different types of sizes: S (small), M (medium) and B (big). HOUSE BRAND PROMOTION-NEW AND DORMANT ACCOUNTS Major Plus Comp Anteriors and Posteriors Super Quality at a super price **Total Retail Value** \$1,932.60 AJOR PLUS COMP Specia Dea 50 Cards Mix & Match BUY **GET** 10 Cards Mix & Match FREE CALL TODAY TO PLACE YOUR ORDER FOR REORDERS **SEE PAGE 20 FOR DISCOUNTS** 805-487-9868

DESCRIPTION

Anterior Set 1x6 MAJOR PLUS COMP

Posterior Set 1x8

MAJOR CHROM

MAJOR PLUS COMP

MAJOR PLUS COMP

WFA Posterior Set 1x8

SHADE GUIDE A1-D4

GROUP/ITEM #

085-T2031

085-T2032

085-T2052

085-S1011





### SPECIALIZED TEETH

# **Major<sup>®</sup> Plus**

### **Anterior Molds and Posterior Molds**

Specifically developed to achieve prestigious functional and aesthetic results in the modern prosthetic. Manufactured with acrylic material, the posterior molds are available in both Anatomical and Technical version. The complete assortment of 48 anterior molds and 26 posterior molds is available in the A1-D4 shade guide.

- High hardness and resistance to abrasion
- Precise occlusal profiles
- New and natural aesthetics of the molds
- Suitable for any mounting technique
- Cuspidal inclination of Anatomical molds between 23° and 25°
- Cuspidal inclination of Technical molds are around 20°
- Cuspidal inclination of High Anatomical molds are 25°
- 4 layers acrylic teeth
- Appropriate for any kind of total and partial prosthesis



GROUP/ITEM #	DESCRIPTION	RETAIL PRICE	SPECIAL PRICE
085-T2011	Major Plus Anterior Set 1x6	\$ 32.21	\$21.45
085-T2012	Major Plus Posterior Set 1x8	\$ 32.21	\$21.45
085-T2026	Major Plus WFA Posterior Set 1x8	\$ 32.21	\$21.45





CALL TODAY TO PLACE YOUR ORDER

805-487-9868

FOR REORDERS

**SEE PAGE 20 FOR DISCOUNTS** 



SEL

Delta∆<sup>™</sup> Anatomic Posterior Teeth

rial & VIVO


These anatomically pleasing 10° posterior teeth are suitable for complete and removable partial dentures. Their carved cusps are very similar in shape to Dentsply Anatoline® teeth. The tooth aesthetics and functionality are complementary with excellent chewing efficiency. This 10° tooth can be used in lingual occlusion configurations or in applications where gentle excursions can offer the patient proper mastication. Imperial Delta $\Delta^{\text{TM}}$  teeth are available in four sizes,  $\Delta$  30X,  $\Delta$  32X,  $\Delta$  34X and  $\Delta$  36X, in all 22 Imperial shades and all 19 Imperial ViVO shades with highly cross-linked intersectional polymer network material with these attributes: maximum wear resistance, vacuum processed, no diatoric retention is necessary, and 100% bondability with methacrylate denture base material.

SPECIAL

\$15.18

GROUP/ITEM #	DESCRIPTION	
083-54-0010	Imperial Delta∆™ Posterior Teeth	\$17.89
083-54-0020	Imperial ViVO <sup>®</sup> Delta∆ <sup>™</sup> Posterior Teeth	\$17.89



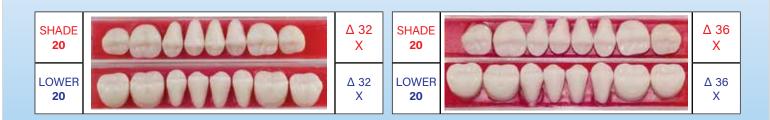

# Available in all 22 Imperial shades and all 19 Imperial ViVO shades

Offer Expires December 31, 2024-not valid with other offers www.americantooth.com | buy online: store.americantooth.com UPPER

Δ34

Х

Δ 34

Х

For Beautiful Smiles



# **Justi Repair Kit Specials**

## Justi® Repair Kit Special

### Justi<sup>®</sup> Teeth for Beautiful Smiles

More Characteristic From the Inventors of Hardened Plastic Teeth Great Discount Offer on Tooth Repair Kits



nerin

For Temporary & Permanent Dentures & Partials

# Justi<sup>®</sup> VP Repair Kit Special

Standard Aesthetic

### A REPAIR KIT that features diverse molds and shades styled after Dentsply's New Hue/Classic tooth lines.

This quality repair kit features a generous selection of 21 cards in an assortment of molds and shades for quick reliable repairs. Recognizing the need for more anterior uppers than lowers, we offer 15 anterior upper cards and 6 anterior lower cards.

Justi VP Anterior Repair Kit 085-56-0014 - \$103.68



**Bleached Shades Repair Kit Special** 



**Justi<sup>®</sup> Teeth<sup>®</sup>** for Beautiful Smiles

From the Inventors of Hardened Plastic Teeth Great Discount Offer on

# Tooth Repair Kits

- 20 cards of Anterior Uppers / Lowers Kit
- 20 cards of Posterior Uppers / Lowers Kit
- 4 popular shades 17, 18, 19, 21
- Restoring a youthful smile

Anterior Kit: 084-50-0302 - \$385.20 Posterior Kit: 084-50-0303 - \$357.80





# Justi<sup>®</sup> Imperial Block Teeth BLOCK TEETH SPECIAL

Made from highly cross-linked blended acrylic material, these block teeth are four posterior teeth joined together on a flat plane. The 0° degree teeth are a time-saving factor in articulation and mounting of the blocks in a denture wax-up. The cuspless carvings eliminate much of the lateral stress during use. Available in Imperial shades 20, 25, 45, 59, 62, 65, 66, 67, 69 and 77, and Bleached shades 17, 18, 19, 21 & Vita Shades. In upper and lower mold 30M only.

GROUP/ITEM	A #	DESCRIPTION		
084-50-0701	1U/Shade #	Justi Imperial Bloc	k Teeth-Uppers (1x8)	\$17.89
084-50-0701	IL/Shade #	Justi Imperial Bloc	k Teeth-Lowers (1x8)	\$17.89
BUY	12	Cards		





Lower left quadrant Lower right quadrant

Now available in 13 of the newest \*Vita type shades. A1, A2, A3, A3.5, A4, B1, B2, B3, C1, C2, C4, D2, D3 and 4 bleached shades

\*Vita Shades' is not a trademark of American Tooth Industries.

# Justi<sup>®</sup> BlockCutter

Hard Copolymer, No Metal!

2 Cards FREE

GET



### COMES IN 22 IMPERIAL<sup>®</sup> AND 7 POPULAR VITA<sup>®</sup> SHADES: A1-A2-A3-A3.5-A4-B2-B3

\*Vita Shades is not a trademark of American Tooth Industries.

GROUP/ITEM #	DESCRIPTION
083-54-0008	Hard copolym

Hard copolymer bar bladed combo 1x16 Comes \$ in one Mold size: 0° - 30M

The process used to control the tooth shading for BlockCutter teeth has been improved. The shade of the bladed mandibulars is now a higher quality match for the maxillary block teeth.

This new block teeth design allows for very fast setups and enhances the power of a removable denture to more efficiently perform mastication with minimal effort. This increased chewing efficiency is done without undue wear on the opposite maxillary plastic teeth.



# Try Justi<sup>®</sup> Mastix<sup>™</sup>

### Hard Copolymer No Metal!



### COMES IN 22 IMPERIAL<sup>®</sup> AND 7 POPULAR VITA<sup>\*</sup> SHADES: A1-A2-A3-A3.5-A4-B2-B3

\*Vita Shades<sup>®</sup> is not a trademark of American Tooth Industries.

GROUP/ITEM #DESCRIPTION083-54-0005Hard copolymer<br/>bladed cusps<br/>combo 1x16

Mastix allows the perfect mastication by the cutting shape of the maxillary bicuspids, the first and second molars, left and right. The cusps are positioned to allow an effective lingualized occlusion, making this arrangement the most powerful for the edentulous patient. This design ensures comfort, functional performance, and aesthetic values without undue wear on the opposing plastic teeth.

mer \$67.90 \$56.76

Offer Expires December 31, 2024–not valid with other offers www.americantooth.com | buy online: store.americantooth.com

### Justi<sup>®</sup> Specialty Teeth Products

For Beautiful Smiles



# **Dr. Levin Lingual Bladed Teeth SPECIAL: BUY 3 GET 1 FREE**

Crafted with aesthetics in mind, the cast metal cross blades that are on the lingual cusps of the maxillary second premolars and molars enable greater mastication. Available in 33° med, 20° med and 33° short in popular BF shades 59, 62, 65, 66, 67, 81, bleached 17, 18, 19, 21 and Vita shades A1, A2, A3, A3.5, A4, B2 and B3.





NOW in Vita<sup>®</sup> Shades and **Bleached Shades** 

**GROUP/ITEM #** 083-54-0001/SHADE 083-54-0003/SHADE

077-65-0057

#### DESCRIPTION

Lingual Bladed Teeth 33°M (1x16) \$136.68 083-54-0002/SHADE Lingual Bladed Teeth 20°M (1x16) \$136.68 Lingual Bladed Teeth 33°S (1x16) \$136.68 Lingual Bladed Tooth Manual\* \$12.59

> \*Also available as a free download at: www.americantooth.com \*Vita Shades\* is not a trademark of American Tooth Industries.

# **Deciduous Starter Kit**

A convenient assortment of Justi Bambino Deciduous teeth, used as spacemaintainers for congenitally missing deciduous teeth and to preserve a space when a primary tooth is exfoliated prematurely. Perfect for labs, orthodontists and pedodontists.



Please specify Anterior or Posterior and Shade. Price: \$270.00

#### SPECIAL Get 10 Bambino<sup>®</sup> Teeth Mix & Match \$10.65 A Card

The Bambino tooth is an aesthetic and functional solution to the problem of the empty spaces left by deciduous teeth.

Made of a copolymer crosslinked resin, this tooth line guarantees impermeability to infiltration of fluids and bacteria, thus preventing discoloration and changing of dimensional character. Studied and developed by orthodontic specialists as a



simplified alternative to orthodontic appliances, these teeth are ideal for use as a space maintainer to ensure a correct line of teeth where both articulation and aesthetic are concerned. Use of these teeth in articular rehabilitations allows for proper mastication, which can only be obtained with perfect articular balance. Available in a singlemold upper and lower anteriors and upper and lower posteriors. Mounted on dark plastic blocks in 2 shades: 59 (milk white), 61 (yellow-straw white).

#### GROUP/ITEM # DESCRIPTION

084-10077*	Justi Bambino Deciduous Upper Anterior Card (1x6)	\$13.50
084-10078*	Justi Bambino Deciduous Lower Anterior Card (1x6)	\$13.50
084-10177*	Justi Bambino Deciduous Upper Posterior Card (1x4)	\$13.50
084-10178*	Justi Bambino Deciduous Lower Posterior Card (1x4)	\$13.50
084-10278*	Justi Bambino Deciduous Combination Card (1x20)	\$54.00
*Please specify	shade	

### Justi<sup>®</sup> Specialty Teeth Products

For Beautiful Smiles

# Justi is your only source in the world for

# Imperial<sup>®</sup> Flatback Plastic Slotted Facings

Justi Imperial Facings are made of the same wear-resistant formula as our Justi Imperial tooth line. The facings are available in a selection of Imperial shades and are pigmented to match the shade guide when used with a standard metallic backing.

The Advantages of Plastic Facings are: 1. Freedom from breakage. Will not chip or fracture while being ground in. Will not break in normal mouth service.

2. Elimination of glazing after grinding to tissue adaptation. Polishing restores the surface so that it is well tolerated by mucosa. Porcelain must be reglazed to prevent irritation.  Adaptability. They may be easily ground to any shape, recarved on the labial or characterized. In case of extreme resorption, pink material may be added to the gums with a pleasing effect.
 Superior Aesthetic. Distinctive anatomy and fidelity of shade. Justi Imperial Facings may be ground to any shape and can be characterized or recarved on the labial to match adjoining teeth. A preformed plastic backing is available and can be easily modified.

Justi<sup>®</sup> Hardened Plastic Teeth comply with ANSI/ADA Specification #15 Published in 2008, identical adoption of ISO Standard #22112:2017

Available in small, medium and large anterior tooth sizes



Single Teeth Available

Ask for a life-like Mold Chart or download at www.AmericanTooth.com

Justi Imperial Flatback Plastic Slotted Facings						
GROUP/ITEM #	DESCRIPTION		PROMO PRICE	PROMO QUANTITY		
084-CL52-0102	Centrals and Laterals (1x4)	\$23.96	\$22.60	When You Buy 5 Mix & Match		
084-CUSP520202	Paired Cuspids (1x4)	\$23.96	\$22.60	When You Buy 5 Mix & Match		
084-52-0104	Facings Shade Guide in 8 Shades	\$13.03	\$11.73	When You Buy 10 pcs.		
084-52-0101	Single Facing (1x1)	\$10.75	\$9.50	When You Buy 5 Mix & Match		
077-52-0105	Facings Backings (50 pieces)	\$49.42	\$36.23	When You Buy 2 Packs		
Special Facing Asso	ortment Kit					
GROUP/ITEM #	DESCRIPTION		PROMO PRICE	PROMO QUANTITY		
084-52-0103	Special Facing Assortment Kit	\$886.32	\$626.75	When You Buy 1 Kit		

# **DISCOUNTS FOR RETAIL CUSTOMERS**



See pages 6-16 for stocking specials. Call your Tooth Express representative for pricing.

Dear Retail Customers,

In order to increase your profits and receive an excellent discount, we are confirming the following price breaks (by tooth lines) for the following lab tooth prices:

MINIMUM SETS (MIX & MATCH) WITHIN TOOTH LINES YOU NEED TO BUY TO RECEIVE DISCOUNT BELOW

Grp #	Description	20-100	101-200	201-300	Over 300
079	Justi Blend Anterior & Posterior (Mix & Match)	10%	15%	20%	
082	Dymon-Hue Anterior & Posterior (Mix & Match)	10%	15%	20%	
085	Dymon-Hue V Anterior & Posterior (Mix & Match)	10%	15%	20%	
083	Imperial and Imperial ViVO Anterior & Posterior (Mix & Match)	10%	15%	20%	Please
084	Bambino Anterior & Posterior (Mix & Match)	10%	15%	20%	ask for a
084	Justi Block Posterior	10%	15%	20%	quotation
085	Major Plus & Major Plus Comp Ant & Post (Mix & Match)	10%	15%	20%	invoice
085	Super Lux Anterior	10%	15%	20%	
085	Major Dent Anterior	10%	15%	20%	
082	Triplex Posterior	10%	15%	20%	

Discounts expire 12/31/24 · No other discounts apply



# Justi Tooth Trays and Cabinets are an ideal storage and inventory system for Justi Teeth. Cabinets are available in three sizes.

- Large will hold up to 4,576 combination Justi anterior and/or posterior tooth cards
- Medium will hold 1,512 anterior and/or posterior cards
- Small has a capacity of 210 to 360 cards of anteriors and/or posteriors

The large counter cabinet does not require inserts. Individual tooth drawer trays contain inventory strips to write the card brand, shade and mold information on, which simplifies reorders. Justi's tray inserts also fit a standard Dentsply 16 1/2'' x 16 7/16'' metal tray, which minimizes brand conversion costs.





### **44-Drawer Tooth Counter**

	Tooth Capacity	Total Cards
Tooth Line	Cards per Drawer per Counter	
Dymon-Hue HPT Anteriors	72	3,168
Dymon-Hue HPT Posteriors	104	4,576
Jusit Imperial Anteriors	80	3,520
Justi Imperial Posteriors	104	4,576
Justi Blend Anteriors	64	2,816
Justi Blend Posteriors	104	4,576

Counter Dimensions: 19" x 19" x 47 1/4"\*

\*Allow an additional 2 inches when casters are installed.

Does not include inserts



with purchase of 200 mix and match Justi Teeth

### **27-Drawer Lab Tooth Counter**

Laminated material; 27 drawers with a capacity of 1,512 anterior and/or posterior cards. Can be used for any Justi tooth brand.

Counter Dimensions: 19" x 19" x 33 1/2"\*

\*Allow an additional 2 inches when casters are installed.

GROUP/ITEM #	DESCRIPTION		SPECIAL
086-55-0201	27-Drawer Lab Tooth Counter	\$2,215.13	
086-55-0201/T	Tooth Counter Trays 10 & 4 Slots	\$33.83	2022.00



### Justi<sup>®</sup> Tooth Cabinets and Tray Specials



### **Lab Tooth Cabinet**

Made of heavy plastic. Ten drawers with 1 tooth insert each. Will hold 320 to 560 anterior and/or posterior tooth cards. Can be used for all Justi tooth brands.

Cabinet Dimensions: 11 1/4" x 13 3/4" x 12 1/2"

 GROUP/ITEM #
 DESCRIPTION

 086-T8030
 Justi Lab Tooth Cabinet 3000
 \$335.43



Since 1864

with purchase of 50 mix and match Justi Teeth

### Justi Lab Tooth Cabinet Drawer Inserts


1	_	-		
-	-	-	-	4
1				
-		-	-	
		1000	_	
in the second				
		1000		
-		-		

**Tooth Drawer Insert** 

Universal Tray Insert Fits 7 anteriors or 7 posteriors 007-65-0015



Tooth Drawer Thick Insert For anterior teeth 007-65-0027

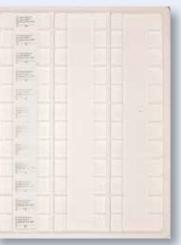

#### Tooth Drawer Thick Insert For posterior teeth 007-65-0028

and the second se		
CONT.		
Converse of		
and the second s		
	800 C	
0.00.00		
Contraction of the second		
100.000		
and the second second		
and the second se		
and the second se		
parts.		
Contract of the second		
and the second second		
「東京市	darden i	
" Andreaded "		
Constant -		
and the second sec		
1		
and the second		
and the second se		
	dimension of the second s	
and the second		
A Strengton		
1 m m m m m m m m m m m m m m m m m m m		
BEE.		
- 10 mm		
	and the second se	
Langer of		
and the second s		
Station .		
and the second s		
and the second se	and the second se	
- BARRIER -		
Safety at		
1		
	 and the second se	
And then		
- Company and		
	and the second se	
1000-		

Drawer inserts shown with optional Tooth Strips.

GROUP/ITEM #	DESCRIPTION	
077-65-0007*	Tooth Drawer Insert - holds 20 sets of the following: Justi Blend Anteriors, ViForm Anteriors and/or Posteriors, Major Plus Comp Posteriors and Triplex Posteriors.	\$ 3.91
077-65-0007/S*	Tooth Strip Insert - for 20-Set Tooth Drawer, anteriors and/ or posteriors.	\$ .75
077-65-0015	Universal Tray Insert - fits 7 Anteriors and/or Posteriors. Fits all Justi tooth cards	\$ 1.64
077-65-0016	Tooth Strips Universal -7 spaces for tooth card information	\$ .66
077-65-0027	Tooth Drawer Thick Insert - for anteriors. Holds 36 sets of the following: Justi Imperial Anteriors, Dymon-Hue HPT Anteriors, Super Lux Anteriors, Bambino and Justi Kind Anteriors and Posteriors, and Major Dent Anteriors.	\$ 3.91

GROUP/ITEM #	DESCRIPTION	
077-65-0027/S*	Tooth Strip Insert - for 36-Set Tooth Drawer, anteriors.	\$ .75
077-65-0028	Tooth Drawer Thick Insert - for posteriors. Holds 28 sets of the following: Justi Imperial Posteriors, Justi Blend Posteriors, Dymon-Hue HPT Posteriors and Justi Imperial Block Posteriors.	\$ 3.91
077-65-0028/S*	Tooth Strip Insert - for 28-Set Tooth Drawer, posteriors.	\$ .75
*Taath Ctuin Incouto	and a shell a subscription.	

\*Tooth Strip Inserts are sold separately.





### **Justi® Denture Base**

- · Short setting
- Long working time
- High color fidelity
- Prevents crazing from solvents



### Justi<sup>®</sup> has the Largest Shade Selection Available with the Most Fibers



BUY 1 25 LBS. DENTURE RESIN POWDER AND GET 1 GALLON OF LIQUID AT NO CHARGE! SAVE \$89.36!

#### SPECIAL - \$621.87

GROUP/ITEM #	DESCRIPTION					
Refills - 25 lb. powd	er					
072-04-0004	Regular pink powder	\$740.44				
072-04-0008	Dark powder	\$740.44				
072-04-0011	B Shade powder	\$740.44				
072-04-0014	072-04-0014 Regular fibered powder					
072-04-0017	Dark fibered powder	\$740.44				
072-04-0020	Clear powder	\$740.44				
072-04-0036	Ethnic Natural	\$740.44				
072-04-0040	Ethnic "A"	\$740.44				
072-04-0044	Ethnic "B	\$740.44				
072-04-0048	Regular orig. fibered	\$740.44				

### Acrylic Materials for prosthetic use

PECIAL - \$61 59



### Justi<sup>®</sup> Denture Base Materials High Impact for Traditional and Microwave Technique

A high-grade methacrylate polymer combined with a cross-linked, solvent-resistant monomer has advantages not found in ordinary denture base materials. It has a short setting and long working time. It is hard, dense and structurally strong, having good model adaptation and high color fidelity. The high-quality cross-linked monomer prevents crazing from solvents such as ether, alcohol, monomer and acetone, etc. This material will bond with any of Justi's tooth lines.

	SPECIAL - Ş	01.55				
GROUP/ITEM #	DESCRIPTION					
Justi Denture Base Kit (20 units)* 1 lb. powder and 8 oz. clear liquid						
072-04-0001	Regular pink powder	\$76.82				
072-04-0006	Dark powder	\$76.82				
072-04-0009	B Shade powder	\$76.82				
072-04-0012	Regular fibered powder	\$76.82				
072-04-0015	Dark fibered powder	\$76.82				
072-04-0018	Clear powder	\$76.82				
072-04-0033	Ethnic Natural	\$76.82				
072-04-0037	Ethnic "A"	\$76.82				
072-04-0041	Ethnic "B"	\$76.82				
072-04-0045	Regular orig. fibered	\$76.82				

	SPECIAL - Ş	45.89
GROUP/ITEM #	DESCRIPTION	
Refills - 1 lb. powde	er	
072-04-0002	Regular pink powder	\$61.12
072-04-0007	Dark powder	\$61.12
072-04-0010	B Shade powder	\$61.12
072-04-0013	Regular fibered powder	\$61.12
072-04-0016	Dark fibered powder	\$61.12
072-04-0019	Clear powder	\$61.12
072-04-0034	Ethnic Natural	\$61.12
072-04-0038	Ethnic "A"	\$61.12
072-04-0042	Ethnic "B"	\$61.12
072-04-0046	Regular orig. fibered	\$61.12

\*Enough material for 20 dentures

Offer Expires December 31, 2024–not valid with other offers www.americantooth.com | buy online: store.americantooth.com

	SPECIAL - \$24	46.33
GROUP/ITEM #	DESCRIPTION	
Refills - 7 lb. bulk po	owder	
072-04-0025	Regular pink powder	\$256.49
072-04-0026	Dark powder	\$256.49
072-04-0027	B Shade powder	\$256.49
072-04-0028	Regular fibered powder	\$256.49
072-04-0029	Dark fibered powder	\$256.49
072-04-0030	Clear powder	\$256.49
072-04-0035	Ethnic Natural	\$256.49
072-04-0039	Ethnic "A"	\$256.49
072-04-0043	Ethnic "B"	\$256.49
072-04-0047	Regular orig. fibered	\$256.49

GROUP/ITEM #	DESCRIPTION	
Refills - Liquid		
072-04-0003	Microlink Liquid 8 oz.	\$15.55
072-04-0031	Microlink Liquid 1 qt.	\$44.70
072-04-0005	Microlink Liquid 1 gal.	\$89.36
072-04-0021**	Micro-Liquid 8 oz.	\$18.80
072-04-0032**	Micro-Liquid 1 qt.	\$62.05
072-04-0022**	Micro-Liquid 1 gal.	\$155.74

\*\* For faster curing use Justi Micro-Liquid

	SPECIAL - \$6	21.87
GROUP/ITEM #	DESCRIPTION	
Refills - 25 lb. powd	er	
072-04-0004	Regular pink powder	\$740.44
072-04-0008	Dark powder	\$740.44
072-04-0011	B Shade powder	\$740.44
072-04-0014	Regular fibered powder	\$740.44
072-04-0017	Dark fibered powder	\$740.44
072-04-0020	Clear powder	\$740.44
072-04-0036	Ethnic Natural	\$740.44
072-04-0040	Ethnic "A"	\$740.44
072-04-0044	Ethnic "B	\$740.44
072-04-0048	Regular orig. fibered	\$740.44

SPECIAL - Free with minimum 7lb. purchase							
GROUP/ITEM # DESCRIPTION							
072-04-0102	Justi Denture Base Shade Guide - 10 Shades	\$14.07					



call today to place your order 805-487-9868



SPECIAL PRICE

### Justi® Pourable Denture Base Material

High performance polymer in a pourable acrylic.

With Justi's Pourable Denture Base Material, you can enjoy a pourable denture base that mixes quickly to produce the perfect consistency required for optimum flow and exceptional performance.

#### Features:

- 4-5 minute pour time and a 20 minute short curing time
- Easy mixing for perfect consistency
- Super-fine particle size enhances strength and tooth bond
- Low porosity and minimal shrinkage
- Cadmium free
- Exceptional color stability



			SPEC	AL PRICE
	GROUP/ITEM #	DESCRIPTION		
	072-50-0001	Pourable db dk fib kit	\$107.78	\$82.00
	072- 50-0002	Pourable db dk fib 1-lb	\$59.70	\$39.80
	072- 50-0003	Pourable db dk fib 7-lb	\$353.08	\$235.88
	072- 50-0004	Pourable db dk fib 25lb	\$1,237.88	\$825.25
	072- 50-0005	Pourable db regular fibered kit	\$107.78	\$82.00
	072- 50-0006	Pourable db regular fibered 1-lb	\$59.70	\$39.80
	072- 50-0007 Pourable db regular fibered 7-lb \$353.0		\$353.08	\$235.38
	072- 50-0008	Pourable db regular fibered 25lb	\$1,237.88	\$825.25
	072- 50-0009	Pourable db original fibered kit	\$107.78	\$82.00
	072- 50-0010	Pourable db original fibered 1lb	\$59.70	\$39.80
	072- 50-0011	Pourable db original fibered 7lb	\$353.08	\$235.38
	072- 50-0012	Pourable db original fibered 25lb	\$1,237.88	\$825.25
	072- 50-0013	Pourable denture base liquid 8oz	\$26.38	\$17.58
	072- 50-0014	Pourable denture base liquid 1qt	\$81.28	\$54.18
	072- 50-0015	Pourable denture base liquid 1gl	\$273.23	\$182.15

### Justi® Repair Material

A universal-type repair material, applicable for repairing all types of acrylic dentures, temporary prostheses and indirect relines. This self-curing methyl methacrylate material is cross-linked and internally pigmented to give a dense, even texture. Justi Repair Material is available in 3 shades of powder: Regular Pink, Regular Fibered Pink and Clear.

Buy 25 lbs.

Get 1 gal. FREE

Buy 7 lbs.

Get 1 gt. FREE

Repair Re	Reparts Sates

		SPECIA	L PRICE				SPECIAL PRICE
GROUP/ITEM #	DESCRIPTION			GROUP/ITEM #	DESCRIPTION		
luoti Bonoir Motor	Justi Repair Material Kit (20 units) 1 lb. powder and 8 oz. liquid				Liquid 8 oz.	\$17.70	\$13.46
Justi nepair iviater	iai Kit (20 units) 110. pov		oz. liquiu	072-06-0015	Liquid 1 qt.	\$46.0	1 \$40.37
072-06-0001	Regular pink	\$96.49	\$65.91	072-06-0013	Liquid 1 gal.	\$105.88	8 \$83.95
072-06-0004	Regular fibered pink	\$96.49	\$65.91	Shade Guide			
072-06-0006	Clear	\$96.49	\$65.91		Repair Material		Free with any
Refills - 1 lb. powd	er			072-06-0100	Shade Guide	\$6.1	2 purchase
072-06-0002	Regular pink	\$78.00	\$52.46				
072-06-0005	Regular fibered pink	\$78.00	\$52.46				
072-06-0007	Clear	\$78.00	\$52.46				
Refills - 7 lb. bulk	powder			Refills - 25 lb. bull	< powder		
072-06-0014	Regular pink	\$331.55	\$242.16	072-06-0010	Regular fibered pink	<\$10	81.71 \$792.06
072-06-0016	Regular fibered pink	\$331.55	\$242.16	072-06-0012	Clear	\$10	81.71 \$792.06
072-06-0017	Clear	\$331.55	\$242.16				

### Justi<sup>®</sup> Tray Material

A self-curing-type acrylic for construction of precision custom trays and stabilized occlusal bite rims. Justi Tray Material is a coarse-type powder that has a high degree of elasticity with a moderate curing time to allow ample working time. After the initial dough stage, the material can be manipulated without sticking to the hands. Having a dead pack, it can be adapted to the model without any significant shrinkage after the cure. The surface roughness aids in the retention of the impression material.

-			SPEC	IAL PRICE
Day Mont	GROUP/ITEM #	DESCRIPTION		
	072-03-1030	30-Unit Professional Package	\$80.30	\$70.04
<b>1</b>	072-03-1025	Bulk Shipment Contents: Powder - 1 25 lb. drum, Liquid - 1 gallon	\$978.56	\$676.20
	Refills:			
10	072-03-0499	Powder - 1 25 lb. drum	\$867.43	\$700.35
	072-03-0699	Liquid - 1 gallon	\$111.11	\$72.45



# Justi<sup>®</sup> IMPAK Pink-Fibered (for Hybrid Dentures)

**Hybrid Denture** – Traditional Hard Acrylic for Denture Teeth Matrix, IMPAK-PF for Soft, Flexible Flanges.

**Stabilized Lower Full Dentures** – Engage previously unusable posterior lingual undercuts.

Pink-Fibered – Matches popular hard Acrylic Shades.
Excellent Bond – Adheres well to Traditional Acrylics.
Adjustable Softness – Alterable Powder/Liquid Ratio.
Added Comfort – Minimal Patient Adjustment period.
Economical Start-Up – No special equipment required.

### PURCHASE OF LIQUID CAN ONLY BE MADE WITH THE POWDER.

GROUP/ITEM #	DESCRIPTION	Intro Price
072-3311	IMPAK Fiber Pink 1lb (454G) Powder Only	\$86.20
072-3312	IMPAK Fiber Pink 6LBS (2.72kg) Powder Only	\$375.80
072-3348	IMPAK Transparent Liquid 11oz (325 CC) Liquid Only Package	\$32.10
072-3306	IMPAK Transparent Liquid 1 Quart (946 CC) Liquid Only Quart	\$51.90



### **Justi® Frame Base Material**

Justi Frame Base Material is a specially formulated microwave-type denture base material used in processing metal framework partial dentures. Can be rapidly cured in a microwave oven. This material can withstand higher heat curing than conventional acrylics. Used only for microwave curing with the Justi VIPI Microflask System. Available in two shades: regular fibered and dark fibered.



GROUP/ITEM #	GROUP/ITEM # DESCRIPTION		SPECIAL PRICE
072-07-0001	Frame Base Regular Fibered -1 lb. Kit - (1 lb. powder, 8 oz. Micro-liquid)	\$72.16	\$52.49
072-07-0002	Frame Base Dark Fibered - 1 lb. Kit - (1 lb. powder, 8 oz. Micro-liquid)	\$72.16	\$52.49



### Justi<sup>®</sup> Crown & Bridge Formula

An outstanding cross-linked polymethyl methacrylate resin to formulate crowns, bridges and inlays. It is a heat-cure method for permanent cases and long-term provisionals that is available in 17 different dentin shades: 20, 25, 30, 35, 40, 45, 50, 55, 60, 62, 65, 66, 67, 69, 70, 77 and 81. It also comes in two enamels: light incisal and grey incisal and a finishing transparent shade. Complete technique instructions are available on request.

		SPECI	AL PRICE				SPECIAL PRICE
GROUP/ITEM #	DESCRIPTION			GROUP/ITEM #	DESCRIPTION		
070-40-0601	1 oz. "SR" Monomer (Solvent Resistant)	\$6.56	\$6.44	070-40-0404/SHADE*	4 oz. Powder - all shades	\$31.05	\$30.1
070-40-0604	4 oz. "SR" Monomer (Solvent Resistant)	\$15.57	\$15.10	070-40-0420/SHADE*	1 lb. Powder - all shades	\$103.64	\$100.5
070-40-0632	32 oz. "SR" Monomer (Solvent Resistant)	\$59.52	\$57.30	070-40-0401/SHADE*	1 oz. Powder - all shades	\$15.57	\$15.1
070-40-0699	1 gal. "SR" Monomer (Solvent Resistant)	\$130.20	\$126.30	070-42-0421	Crown and Bridge Formula (SR Shade Selector)	\$25.86	Free with any

\*Specify Shade





### **Justi® Imperial Namilon**

An internally pigmented, cross-linked methyl methacrylate resin for permanent crown and bridge restorations. It is ideal for quality, highstrength, long-term provisional cases such as implants and full-mouth reconstructions. Namilon is to be used with the crown and bridge investment procedures, i.e., wax-up, investment, boil-out and packing. Because Namilon is internally pigmented, it will not streak, bleed or change color. Namilon features an extended working time, and once polymerized in heat, it cures quickly. Restorations are durable, solvent and wear resistant.

		SPECI	AL PRICE
GROUP/ITEM #	DESCRIPTION		
070-20-1000	Namilon Kit - Contents:16 - 1 oz. bottles Namilon Powder in Imperial Shades: 60, 62, 64, 65, 66, 67, 68, 69, 70, 77, 81, 82, Incisal A, Incisal B, Opaque and Trans- parent, • 2 - 1 oz. bottles Liquid Resin, 2 - Mixing Jars, • 1 - Measuring Spoon, 1 - Glass Pipette, • 1 - Imperial Shade Guide, • 1 - Namilon Instructions	\$378.37	\$367.02
070-20-0401/SHADE*	1 oz. Powder - all shades	\$21.46	\$20.83
070-20-0404/SHADE*	4 oz. Powder - all shades	\$32.43	\$31.46
070-20-0420/SHADE*	16 oz. Powder - all shades	\$111.17	\$107.84
070-20-0601	Namilion 1 oz. Liquid	\$8.38	\$8.14
070-20-0604	Namilion 4 oz. Liquid	\$16.33	\$15.85
070-20-0608	Namilion 8 oz. Liquid	\$27.75	\$26.93
070-20-0632	Namilion 32 oz. Liquid	\$61.40	\$60.23

\*Specify Shade



### Justi<sup>®</sup> Soft Denture Reliner Semi-Rigid, Soft and Very Soft, Heat Cure

An ethyl methacrylate-type material that can be heat cured traditionally or processed in the microwave oven within minutes when using the Justi Vipi Microflask system. This pink soft liner blends well with most denture base materials. It has a choice of three textures- semi-rigid, soft and very soft. The material is tissue tolerable and will improve mouth comfort from spiny ridges or other causes of chronic soreness. Doctors, technicians and patients will all be exceptionally pleased with this high-quality and long-lasting material.

		SPEC	IAL PRICE
GROUP/ITEM #	DESCRIPTION		
072-14-0001	Justi Denture Reliner Kit - Semi-rigid Contents 1 - 8 oz. Liquid - 1 lb. powder	\$91.10	\$78.49
072-14-0002	Justi Denture Reliner Kit - Soft Contents 1 - 8 oz. Liquid - 1 lb. powder	\$91.10	\$78.49
072-14-0003	Justi Denture Reliner Kit - Very soft Contents 1 - 8 oz. Liquid - 1 lb. powder	\$91.10	\$78.49
072-14-0004	Justi Denture Reliner - 8 oz. Liquid - Semi-rigid	\$30.45	\$25.36
072-14-0005	Justi Denture Reliner - 8 oz. Liquid - Soft	\$30.45	\$25.36
072-14-0006	Justi Denture Reliner - 8 oz. Liquid - Very Soft	\$30.45	\$25.36
072-14-0007	Justi Denture Reliner - 1 lb. Powder	\$60.59	\$43.47



Item # 070-28-0020

### Justi<sup>®</sup> Pakto System "The Artesians and Ceramist's Friend"

Same technique as is used with porcelains. The Justi Pakto® System is a high-quality resin-to-metal alternative to porcelain.

Justi Pakto is a crown and bridge resin material system designed to save time and labor while offering maximum physical properties to simulate the hardness of natural dentition. Justi Pakto is not a methyl methacrylate, a composite or a porcelain. It is 100% poly glycol dimethacrylate resin and is processed in a dry-heat pressurized chamber.

The areas of application are:

- A) Crown and bridge veneers;
- B) Facings on partial dentures and attachment frames;
- C) Inlays of all types;
- D) Maryland-type bridges; and
- E) Temporary crowns and bridges.

Justi Pakto can be used with any ADA-recommended metal with no color change or fractures. The opaque has superior masking qualities and chemically bonds to the tooth-shaded resin material.

Color verification can be done before the resin is cured since the shade is identical when cured. Justi Pakto is shaded to match the Vita-Lumin<sup>®</sup> Shade System. (Vita-Lumin is not a registered trademark of American Tooth Industries.)

The same tools used to process porcelain can be used for Justi Pakto. The Justi Pakto Curing System uses a 4-quart, dry-heat pressurized pot.

Replacement parts sold seperately.

GROUP/ITEM # DESCRIPTION					
Pakto Pressure Pot System:					
CONTENTS: 1 - Pressure Pot - 4 quarts, 070-28-0020 1 - Curing Vessel, 2 - Security Spring Clips, 1 - Air Supply Hose (3-month warranty on system)	\$545.01				
Justi Pakto Kit:					
<ul> <li>16 - 1 oz. Bottles Justi Pakto Powder in Vita* Shades: A1, A2, A3, A4, B2, B3, B4, C1, C2, C3, C4, D2, D3, Body Shades, Transparent Overlay and 2 Bottles Incisal, 8 - 1/4 oz. Bottles Justi Pakto Modifying Stains: white, brown, yellow, gray, blue, pink, orange and dark orange, 2 - 1/4 oz. Bottles of Justi Pakto Liquid, 1 - 1/4 oz. Reten- tion Beads, 1 - 1/4 oz. Bead Adhesive, 2 - Dappen Dishes, 1 - Glass Pipette, 1 - Plastic Scoop, 1 - Instruction Manual. (Dappen dishes, eye dropper and plastic scoop are not available for reorder)</li> <li>1 - Justi Pakto Shade Guide - 13 shades, double-sided blue frame</li> </ul>	\$684.05				
Justi Pakto Refills:					
070-28-0002/SHADE 1 oz. Pakto Powders - All Shades (Includes Body Shades, Incisal or Transparent)	\$30.19				
4 Delite Devidere - All Cherder (Jerebider Devide Cherder	\$108.26				
070-28-0003/SHADE 4 oz. Pakto Powders - All Shades (Includes Body Shades, Incisal or Transparent)					
	\$9.25				

GROUP/ITEM #	DESCRIPTION	
Pakto Pressure Pot	System Replacement Parts:	
072-28-0021	Curing Vessel Pot	\$30.04
072-28-0022	Security Spring Clip	\$26.12
072-28-0023	Air Supply Hose	\$30.20
072-28-0026	Gasket for Pakto Pot	\$41.88
072-28-0027	Replacement Handle for Pakto Pot - Top	\$26.75
072-28-0027/B	Replacement Handle for Pakto Pot - Bottom	\$26.75
Justi Pakto Shade (	Guide:	
070-28-0030	Justi Pakto Shade Guide	\$20.38
070-28-0032	Justi Pakto Shade Guide - 13 shades, double- sided blue frame	\$35.17
Justi Pakto Supplie	s:	
070-28-0006	1 oz. Pakto Liquid	\$20.38
070-28-0007	4 oz. Pakto Liquid	\$72.34
070-28-0008	1/4 oz. Retention Beads	\$10.94
070-28-0009	1/4 oz. Bead Adhesive	\$10.41
070-29-0001	1 oz. Justi Pakto Pure Alloy	\$28.96
070-29-0002	3 oz. Justi Pakto Pure Alloy	\$33.04
070-29-0001K	1 Kilo Justi Pakto Pure Alloy	\$1019.69
077-66-0003	Pakto Technique Video	\$61.31

Vita-Lumin is not a registered trademark of American Tooth Industries.

### Resin for prosthetic use









### Justi<sup>®</sup> Temporary Crown & Bridge Resin

### 3 Shades - Light, Medium and Dark

An autopolymerizing self-cure powder and liquid resin for fabricating a crown, bridge or splint. It is very hard, color stable and easy to polish and comes in Light (shade 59), Medium (shade 62), and Dark (shade 69). This material has a low exothermic heat peak and therefore, is tolerable to the oral cavity; it does not burn.

		SPECI	AL PRICE		
GROUP/ITEM #	DESCRIPTION				
072-25-0001	Professional Kit - Boxed, 3 oz. each Light, Medium and Dark powders and 8 oz. liquid	\$82.77	\$80.29		
Refills:					
072-25-0002	Light Powder (85 grams) 3 oz.	\$19.51	\$18.93		
072-25-0003	Medium Powder (85 grams) 3 oz.	\$19.51	\$18.93		
072-25-0004	Dark Powder (85 grams) 3 oz	\$19.51	\$18.93		
072-25-0005	Clear Liquid (240 ml) 8 oz.	\$24.53	\$23.81		

### Justi<sup>®</sup> Tooth Shades

Cold Cure Temporary Crown and Bridge Acrylic for Labs is a specially prepared formulation for dental laboratories to construct and build strong temporary crowns and bridges. Crosslinking monomer with high craze resistance and toughness. Kit available in Imperial/Bioform type shades 59, 62, 65, 66, 67, 69 and 81. New Bleached shades 17, 18, 19 and 21 are also available in 1lb. bottles.

			SPECIAL PRICE
GROUP/ITEM #	DESCRIPTION		
072-25-0010	Complete Tooth Shade Kit Contents: 7 - 3 oz. Tooth Shade powders, 2 - 8 oz. liquid bottles.	\$250.21	\$172.67
072-25-0011/ SHADE*	Imperial Bleached Shade Powder Contents: 1 lb. bottle of tooth shade powder (Be sure to indicate the shade – example 072-25-0011/18.)	\$113.93	\$78.49
Refills available:			
072-25-0008	8 oz. bottle of liquid	\$17.85	\$12.25
072-25-0009/ SHADE*	3 oz. bottle of tooth shade powder (Include shade number with request.)	\$20.72	\$17.08

\*Specify Shade

### Justi<sup>®</sup> Study Model Material

A methyl methacrylate self-curing material that can be poured into a silicone or rubber mold with a recommended technique that will duplicate the model in an off-white/cream color. It can be used for fabricating study model cases that emphasize the appliance rather than the model. This material can be used for many other applications where an off-white resin is needed with an autopolymerization system. Cure time is approximately 20 minutes in open air. Pressure pot technique is recommended for a more compact, dense model.

GROUP/ITEM #	DESCRIPTION		
072-13-0001	Study Model Kit (1 lb.) 1 lb. powder (450gr.) 8 oz. Liquid (240ml.)	\$69.67	\$48.51
072-13-0002	Study Model Kit (5 lb.) 5 lb. powder (2250 gr.) 1 qt. Liquid (960ml.)	\$279.19	\$185.32
072-13-0003	Study Model Powder 1 lb. powder (450gr.)	\$56.09	\$46.00
072-13-0004	Study Model Liquid 8 oz. liquid (240ml.)	\$13.59	\$11.21







### Justi<sup>®</sup> Opaquer

A resin opaquer for masking metal in veneer crowns with outstanding quality of adhesion. It is also very helpful in holding veneer color. Primarily used with Namilon or Justi Crown and Bridge Formula. However, it can also be used with other polymerbased crown and bridge materials. Brown and pink opaquers are used when experience shows they are needed.

•	5	SPECI	AL PRICE
GROUP/ITEM #	DESCRIPTION		
070-30-1000	Justi Opaquer Kit Contents: 2 - 1/2 oz. Liquid / 5 - 1/4 oz. Powder in shades white, gray, yellow, brown and pink / 1 - Sable Application Brush / 1 - Glass Pipette / 1 - Polyethylene mixing slab / Instructions for use	\$68.77	\$66.72
070-30-0400/SHADE	Powder Refills - All Shades - 1/4 oz.	\$9.76	\$9.48
070-30-0600	Liquid Refills - 1/2 oz.	\$9.10	\$8.83

### Justi<sup>®</sup> Resin Cement -Best Cement on the Market For Permanent Cementation

A unique copolymer offering greater adhesive properties and insolubility. Justi Resin Cement offers a lower film thickness and is radio-opaque.

### **Features:**

- Shows adhesion to tooth structure superior to any other cement tested after 24 hours immersion in water.
- Meets the requirement for film thickness in ADA Specification #8, when used according to instructions.
- Available in 6 shades: Natural Yellow, Natural Gray, White, Pink, Brown and Clear.
- Insoluble in mouth fluids.
- As harmless to the pulp as present resin filling materials.
- Radio-opaque.

		SPECIAL	PRICE
GROUP/ITEM #	DESCRIPTION		
070-80-1000/SHADE	Justi Resin Cement Kit Contents: 2 - 1/2 oz. Powder; 1 - 1/2 oz. Liquid; 1 - Measuring Spoon, 1 - Eye Dropper and Instructions for Use	\$42.34	\$41.08
070-80-0400/SHADE	Powder Refill (1/2 oz.) and a Measuring Spoon Please specify shade when ordering.	\$15.97	\$15.50
070-80-0600	Liquid Refill (1/2 oz.) and Eye Dropper	\$10.38	\$10.07

## Justi<sup>®</sup> Resin Modifying Stains

Justi Resin Modifying Stains are primarily used for modifying Justi Namilon shade powders and Justi Crown and Bridge formula. However, they can be used with other types of methacrylates. Also used for characterization of the cervical areas, for decalcification marks and check-lines in plastic teeth, etc. Available in shades: yellow, gray, white, brown, pink, blue, light orange and dark orange.

		SFECIAL	FNICE
GROUP/ITEM #	DESCRIPTION		
070-60-1000	Justi Modifying Stain Kit - includes 1/4 oz. of each shade	\$99.77	\$71.77
070-60-0400/SHADE	Refills - 1/4 oz. (Please specify shade when or- dering)	\$9.10	\$12.59



### Justi<sup>®</sup> Cold-Cure Ortho Resin

An autopolymerizing material primarily used in the construction of orthodontic appliances. The powder is clear and when used with the transparent liquid, will produce a clear base material. A pink liquid is also available, which when used with the clear powder, will produce a translucent pink base material. The powder is a low-flow polymer that absorbs the monomer rapidly. This prevents the material from slumping when using the liquid dropper technique.

			SPECIAL PRICE
GROUP/ITEM #	DESCRIPTION		
Justi Ortho Resin Kit	- 1 lb. powder, 8 oz. liquid, eye dropper and 3 oz. jet dispenser		
072-05-0001	Ortho Resin Kit w/Pink Liquid	\$86.96	\$78.49
072-05-0006	Ortho Resin Kit w/Clear Liquid	\$86.96	\$78.49
072-05-0014	Ortho Resin Kit w/Crystal Liquid	\$89.03	\$80.30

		SPEC	CIAL PRICE
GROUP/ITEM #	DESCRIPTION		
LIQUID REFILLS			
072-05-0003	Pink - 8 oz.	\$18.32	\$12.95
072-05-0009	Pink - 1 qt.	\$51.54	\$39.20
072-05-0005	Pink - 1 gal.	\$109.03	\$78.71
072-05-0007	Clear - 8 oz.	\$18.32	\$12.95
072-05-0010	Clear - 1 qt.	\$51.54	\$39.30
072-05-0008	Clear - 1 gal.	\$109.03	\$78.71
072-05-0015	Crystal - 8 oz.	\$22.30	\$15.35
072-05-0016	Crystal - 1 qt.	\$53.08	\$41.71
072-05-0017	Crystal - 1 gal.	\$121.59	\$83.71

			SPECIAL PRICE
GROUP/ITEM #	DESCRIPTION		
POWDER REFILLS			
072-05-0002	1 lb.	\$67.91	\$52.08
072-05-0002/7	7 lb.	\$311.79	\$251.26
072-05-0004	25 lb.	\$962.84	\$742.28





### Pozzi<sup>®</sup> Spruzzo™

### Colored Acrylic Powders - Solid and Glitter Colors

Pozzi blended colored powders were designed for use with our clear monomer. These colors were created for those wanting a more solid color appearance. Designed for use with the salt-and-pepper technique, these colors are sure to catch every teenager's eye! Each color comes in a 4-ounce easy-to-use spray bottle and is available in many popular colors.

			SPEC	CIAL PRICE
GROUP/ITEM #	DESCRIPTION			
040-40-0034	Spruzzo Powder 4 oz.	Yellow	\$16.39	
040-40-0035	Spruzzo Powder 4 oz.	White	\$16.39	
040-40-0036	Spruzzo Powder 4 oz.	Orange	\$16.39	
040-40-0037	Spruzzo Powder 4 oz.	Black	\$16.39	\$12.42
040-40-0038	Spruzzo Powder 4 oz.	Brown	\$16.39	
040-40-0039	Spruzzo Powder 4 oz.	Blue	\$16.39	
040-40-0040	Spruzzo Powder 4 oz.	Red	\$16.39	
040-40-0040/SG	Spruzzo Shade Guide		\$31.07	\$21.75
040-40-0050	Spruzzo Glitter Powder 4 oz.	Yellow Glitter	\$34.39	
040-40-0051	Spruzzo Glitter Powder 4 oz.	Green Glitter	\$34.39	¢05.05
040-40-0052	Spruzzo Glitter Powder 4 oz.	Clear Glitter	\$34.39	\$25.95
040-40-0053	Spruzzo Glitter Powder 4 oz.	Hot Pink Glitter	\$34.39	
				ψΞοι.

### Pozzi<sup>®</sup> Ortho Color Liquids

Use of these colors results in beautiful and unique appliances. Colors are dense and bright and are manufactured by our factory in Oxnard, California.



### **Color Guides**

A display of 13 fun colors to assist patients in choosing their color.



			SPECIAL -	\$7.25
	GROUP/ITEM #	DESCRIPTION		
	041-RP6031-01	1 oz. Pozzi Ortho Color Liquid "C"		\$9.23
	041-RP6122-01	1 oz. Pozzi Ortho Color Liquid "A"	Clear	\$9.23
	041-RP6123-01	1 oz. Pozzi Ortho Color Liquid "B"	Trans. Pink	\$9.23
ш	041-RP6124-01	1 oz. Pozzi Ortho Color Liquid "D"	Red	\$9.23
	041-RP6125-01	1 oz. Pozzi Ortho Color Liquid "E"	Purple	\$9.23
	041-RP6126-01	1 oz. Pozzi Ortho Color Liquid "F"	Magenta	\$9.23
	041-RP6127-01	1 oz. Pozzi Ortho Color Liquid "G"	Orange	\$9.23
OUNCE	041-RP6128-01	1 oz. Pozzi Ortho Color Liquid "H"	Fluor. Yellow	\$9.23
-	041-RP6129-01	1 oz. Pozzi Ortho Color Liquid "I"	Green	\$9.23
	041-RP6130-01	1 oz. Pozzi Ortho Color Liquid "L"	Turquoise	\$9.23
	041-RP6131-01	1 oz. Pozzi Ortho Color Liquid "M"	Blue	\$9.23
	041-RP6132-01	1 oz. Pozzi Ortho Color Liquid	Black	\$9.23
	041-RP6133-01	1 oz. Pozzi Ortho Color Liquid	Gold	\$9.23

			SI LOIAL - Ş	24.15
	GROUP/ITEM #	DESCRIPTION		
	041-RP6031-04	4 oz. Pozzi Ortho Color Liquid "C"		\$29.48
	041-RP6122-04	4 oz. Pozzi Ortho Color Liquid "A"	Clear	\$29.48
	041-RP6123-04	4 oz. Pozzi Ortho Color Liquid "B"	Trans. Pink	\$29.48
щ	041-RP6134-04	4 oz. Pozzi Ortho Color Liquid "D"	Red	\$29.48
вотт	041-RP6135-04	4 oz. Pozzi Ortho Color Liquid "E"	Purple	\$29.48
	041-RP6136-04	4 oz. Pozzi Ortho Color Liquid "F"	Magenta	\$29.48
ЦЩ Ц	041-RP6137-04	4 oz. Pozzi Ortho Color Liquid "G"	Orange	\$29.48
OUNCE	041-RP6138-04	4 oz. Pozzi Ortho Color Liquid "H"	Fluor. Yellow	\$29.48
4	041-RP6139-04	4 oz. Pozzi Ortho Color Liquid "I"	Green	\$29.48
	041-RP6130-04	4 oz. Pozzi Ortho Color Liquid "L"	Turquoise	\$29.48
	041-RP6141-04	4 oz. Pozzi Ortho Color Liquid "M"	Blue	\$29.48
	041-RP6142-04	4 oz. Pozzi Ortho Color Liquid	Black	\$29.48
	041-RP6143-04	4 oz. Pozzi Ortho Color Liquid	Gold	\$29.48

			SPECIAL -	\$29.28
	GROUP/ITEM #	DESCRIPTION		
	041-RP6031-08	8 oz. Pozzi Ortho Color Liquid "C"		\$49.41
	041-RP6122-08	8 oz. Pozzi Ortho Color Liquid "A"	Clear	\$49.41
	041-RP6123-08	8 oz. Pozzi Ortho Color Liquid "B"	Trans. Pink	\$49.41
щ	041-RP6124-08	8 oz. Pozzi Ortho Color Liquid "D"	Red	\$49.41
BOTTLE	041-RP6125-08	8 oz. Pozzi Ortho Color Liquid "E"	Purple	\$49.41
	041-RP6126-08	8 oz. Pozzi Ortho Color Liquid "F"	Magenta	\$49.41
OUNCE	041-RP6127-08	8 oz. Pozzi Ortho Color Liquid "G"	Orange	\$49.41
lo lo	041-RP6128-08	8 oz. Pozzi Ortho Color Liquid "H"	Fluor. Yellow	\$49.41
œ	041-RP6129-08	8 oz. Pozzi Ortho Color Liquid "I"	Green	\$49.41
	041-RP6130-08	8 oz. Pozzi Ortho Color Liquid "I"	Turquoise	\$49.41
	041-RP6131-08	8 oz. Pozzi Ortho Color Liquid "M"	Blue	\$49.41
	041-RP6132-08	8 oz. Pozzi Ortho Color Liquid	Black	\$49.41
	041-BP6133-08	8 oz. Pozzi Ortho Color Liquid	Gold	\$49.41

			SPECIAL - \$	68.14
	GROUP/ITEM #	DESCRIPTION		
	041-RP6031-16	1 pint Pozzi Ortho Color Liquid"C"		\$98.67
	041-RP6122-16	1 pint Pozzi Ortho Color Liquid"A"	Clear	\$98.67
	041-RP6123-16	1 pint Pozzi Ortho Color Liquid"B"	Trans. Pink	\$98.67
	041-RP6134-16	1 pint Pozzi Ortho Color Liquid"D"	Red	\$98.67
BOTTLE	041-RP6135-16	1 pint Pozzi Ortho Color Liquid"E"	Purple	\$98.67
ЮТ	041-RP6136-16	1 pint Pozzi Ortho Color Liquid"F"	Magenta	\$98.67
PINT E	041-RP6137-16	1 pint Pozzi Ortho Color Liquid"G"	Orange	\$98.67
Ы	041-RP6138-16	1 pint Pozzi Ortho Color Liquid"H"	Fluor. Yellow	\$98.67
	041-RP6139-16	1 pint Pozzi Ortho Color Liquid"I"	Green	\$98.67
	041-RP6130-16	1 pint Pozzi Ortho Color Liquid"I"	Turquoise	\$98.67
	041-RP6141-16	1 pint Pozzi Ortho Color Liquid"M"	Blue	\$98.67
	041-RP6142-16	1 pint Pozzi Ortho Color Liquid	Black	\$98.67
	041-RP6143-16	1 pint Pozzi Ortho Color Liquid		\$98.67

SPECIAL - \$24.15

### Accessories for denture maintenance





DENTURE CLEANER

Renew is a Premium Denture Cleaner in a class of its own. Renew powder is a concentrated denture cleaner with a formulation much more powerful than effervescent denture tablets and powders currently on the market. This formulation allows for superior cleaning in a shorter period of time. Renew has professional strength denture cleaning and is backed by a money back guarantee. It cleans better than others, removing: tobacco, iron, plaque, coffee, calcium, calculus, tea, tarter and bacteria. Renew is sold in a generous 200 gram bottle, which provides 50 denture cleanings.

Each case has 50 free samples, making it easier for patients to try Renew free of charge.

Item# 6261030015Renew Denture Cleaner Bottle 200 grams\$12.37Item# 6261030022Renew Denture Cleaner Sample Pack\$0.63



# **Renew Display Kit**

Renew has been on the market for 20 years in Canada and is only sold by dental professionals. Each case has 50 free samples, making it easier for patients to try Renew free of charge. Since Renew is not available in retail stores, having the bottles displayed and readily available will encourage purchase right at your office.

### S495.08 GROUP/ITEM # 6261030 KIT 6261030 BOX Re Celoso BOX Re

**Kit Special** 

# UP/ITEM #DESCRIPTION030 KITRenew Denture Cleaner Kit/48030 BOXRenew Display Box\$18.52

BUY

GE1

# JUSTI DENTURE BRUSH

Just like natural teeth, for optimal denture hygiene, denture brushing is essential. Daily brushing with the Justi Denture Brush is the safest, most effective way to keep a denture looking its best. Proper brushing makes a denture last longer and helps to keep the mouth clean and fresh. The design of the brush enables complete removal of deposits and food residues.

### Essential Denture Cleaning...

The rounded nylon bristle tips will not harm denture materials. The large ergonomic handle provides positive brush movement control for easy denture cleaning. For more thorough cleaning, combine brushing with periodic soaking in Renew Denture Cleaner.

GROUP/ITEM # DESCRIPTION 077-66-0004 1 Justi Denture Brush - Blue



1 Case 48-Bottle


**1 Display Free** 



# **Removable Prosthodontics Manual**

#### A "must" for professional dental technicians

This comprehensive manual is for technicians who want to study various fabrication procedures used in the dental lab and the basic science needed to successfully communicate with dentists. It is also recommended for technicians who are preparing to take the certification examinations. They are written by U.S. Air Force Dental Laboratory Instructors.

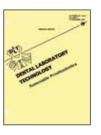

#### GROUP/ITEM # DESCRIPTION 077-65-0031 Removable Prosthodontics - Vol II \$60.31

SPECIAL PRICE \$44.79

# Justi<sup>®</sup> Denture Manual

### Learn How to Fabricate Dentures Using Justi Teeth

A manual for the selection, articulation and fabrication of dentures using Justi plastic teeth. An invaluable guide for the dental laboratory technician. State-of-the-art illustrations describing some of the types of removable prosthodontic appliances and techniques used in dentistry today. An invaluable aid to both teaching and learning.



### AVAILABLE FREE ONLINE

**GROUP/ITEM # DESCRIPTION** 077-65-0006 Justi Tooth Manual \$43.11

# **Denterms Dental Dictionary**

Terminology for the Dental Team-A glossary of dental terminology slated toward those terms associated with prosthetic dentistry. Recommended use by schools of dental laboratory technology and hygiene.

			SPECIAL	5
GROUP/ITEM #	DESCRIPTION		PRICE	fi)
077-65-0049	Denterms- Dental Dictionary	\$35.47	\$11.49	0

For free download go to store.americantooth.com

# Justi<sup>®</sup> Caliper

A measuring device used for measuring parameter landmarks on a working model or facial dimensions to aid in tooth mold selection.

GROUP/ITEM #		SPECIAL PRICE
077-65-0034	\$18.17	\$12.08

Offer Expires December 31, 2024-not valid with other offers www.americantooth.com | buy online: store.americantooth.com



# Justi<sup>®</sup> Papilla Gauge

**Denture Measuring Device** The permanent position of the papilla can provide great diagnostic value in reproducing tooth position in the construction of a denture. A measuring device that uses this landmark will record, with accuracy, the horizontal and vertical plane to define the position of the incisal edges and the central incisors. These measurements can be recorded, along with the mold and shade, on the patient's records for present or future application. For more information, please contact us.



**GROUP/ITEM # DESCRIPTION** 077-65-0013 Justi Papilla Gauge \$125.91

SPECIAL PRICE \$92.00 ea.

# Justi<sup>®</sup> Utility Template

The Justi Template is a quality plastic occlusal plane guide tool specifically designed for the needs of today's productive technicians for tooth setup and occlusal bite rim fabrication. A guide relation of 80 degree angled sides are calibrated with a millimeter rule guide measuring 0 to 70 mm with markings at 18 and 22 mm.

**GROUP/ITEM # DESCRIPTION** 077-65-0041A Justi Utility Template

SPECIAL PRICE \$12.22 \$8 11

# **Technician's Utility Ruler**

The Justi Technician's Utility Ruler will aid the technician in the most common situations: fabrication of occlusal rims, verification of the maxillary central incisor position, determination of normal or crossbite situations, determination of posterior molds, measuring for the six anterior teeth and related information. Comes complete with instructions for use.



#### GROUP/ITEM # DESCRIPTION

077-65-0012 Technician's Utility Ruler \$5.58 Buy 10 for \$43.59

# Justi<sup>®</sup> Pocket Lens

### **Details, Details, Details!**

You'll love how clearly details can be seen with our new Pocket Lens. The applications around your lab or dental office are limitless.

GROUP/ITEM #	DESCRIPTION	
077-65-0035	Justi Pocket Lens	\$4.42



CALL TODAY TO PLACE YOUR ORDER 805-487-9868





Offer Expires December 31, 2024–not valid with other offers www.americantooth.com | buy online: store.americantooth.com

# Alginates for impressions

FDA 510(K) Clearance Approval

High precision alginate for impressions with chromatic phase indicator. Innovative advanced formula for tropical climates and hot environments.

Extremely refined formulation. Thanks to the MAP technology (modified atmosphere packaging) the product shelf life is 6 years.

Chromatic indication of the working phases for easy and precise results, regardless of water temperature and hardness. Setting time: 3'20" (at 23°C).

WORKING/PACKING

SETTING

\$81.87

Excellent gypsum compatibility. Very high compressive strength and recovery from deformation.

recovery from deformation.

• Excellent gypsum compatibility. Very high compressive strength and

Excellent biocompatibility.

#### Setting Times

Mixing Time: 35" Working Time: 2'00" Minimum Setting Time in mouth: 1' 20" Initial Setting Time: 2'00" Total Setting Time: 3'20"

003-A5024-4 Alginplus Tropical Box of 4 lbs.

**ALGINPLUS TROPICAL** 

Superior elasticity and tear resistance. Stabilized hypoallergenic formula. High flow and easy mixing.

Dust free powder for safe working and

CHROMATIC PHASES ALGINPLUS MIXING

## Accessories for impressions

# NYLON SPATULAS FOR ALGINATES

003-2803/BS	Alginmax Blue Spatula	\$2.09
003-2803/WS	Alginmax White Spatula	\$2.09

# **ISOLMAJOR**

Universal plaster/resin separator

L2710 Isolmajor liquid	500 ml	\$24.48
------------------------	--------	---------









**Main Features** 

.

• Aroma: Fruits.

Color: Yellow.

cleaning ease.







# **ORMAPLUS LAB 70**

Vinyl polysiloxane based, addition-curing material, exclusively designed for extra oral use in the dental laboratory.

#### TECHNICAL DATA

- Hardness (after 24h): 70 (±3) Shore-A.
- Recovery from deformation: 99,5%
- Mixing ratio: 1:1
- Linear dimensional change: 0,05%
  Perfect compability with self and heat curing acrylic resins.

High heat resistance without alterations.High precision in detail reproduction.

d heat and bridges. • CAD compatibility.

INTENDED USE

Duplications.

Prosthesis repair.

· Matrices for fixed and total prosthesis.

Matrices for temporary crowns

 M4300
 Ormaplus Lab 70 (0,75 kg Base + 0,75 kg Catalyst)
 \$117.85

 M4310
 Ormaplus Lab 70 (5 kg Base + 5 kg Catalyst)
 \$375.00



# **ORMAPLUS LAB 90**

Vinyl polysiloxane based, addition-curing material, exclusively designed for extra oral use in the dental laboratory.

#### **TECHNICAL DATA**

- Hardness (after 24h): 90 (±3) Shore-A.
- Recovery from deformation: 99,5%
- Mixing ratio: 1:1
- Linear dimensional change: 0,05%
- Perfect compability with self and heat curing acrylic resins.
- High heat resistance without alterations.
- High precision in detail reproduction.

#### INTENDED USE

- Matrices for fixed and total prosthesis.
- Duplications.
- Prosthesis repair.
- Matrices for temporary crowns and bridges.
- CAD compatibility.

M400	00	Ormaplus Lab 90 (0,9 kg Base + 0,9 kg Catalyst)	\$109.79
M40	10	Ormaplus Lab 90 (5 kg Base + 5 kg Catalyst)	\$321.44





#### **Condensation Curing Silicones**





### **ORMALAB 85**

· Matrices for fixed and total prosthesis.

Matrices for temporary crowns and bridges.

A polysiloxane based, condensation curing material, exclusively designed for extra oral use in the dental laboratory. Ormalab 85 is indispensable whenever high quality performances are required.

INTENDED USE:

· Prosthesis repair.

CAD reproduction

- Long lasting dimensional stability.
- Very high compressive strenght and recovery from deformation.

M8501	Ormalab 85 (1600 g / 900 ml)	\$58.88
M8525	Ormalab 85 (5 kg)	\$152.05
M8520	Ormalab 85 (5 kg) + 2 Ormactivator (60 ml)	\$222.38

### **ORMALAB 75**

A polysiloxane based, condensation curing material, exclusively designed for extra oral use in the dental laboratory. Ormalab 75 is useful whenever high quality performances are required.

#### **TECHNICAL DATA**

- Hardness (after 24h): 75(±3) Shore-A.
- Perfect compatibility with cold curing
- acrylic resins.
- High precision in details reproduction.Extremely high heat resistance:
- over 120°C.
- Long lasting dimensional stability.Very high compressive strength and
- recovery from deformation.

#### INTENDED USE

- Matrices for fixed and total prosthesis.
- Prosthesis repair.
- Matrices for temporary crowns and
- bridges.
- CAD compatibility.

M7500	Ormalab 75 (1600 g / 900 ml)	\$64.56
M7510	Ormalab 75 (5 kg)	\$157.73
M7520	Ormalab 75 (5 kg) + 2 Ormactivator (60 ml)	\$231.13



### **ORMALAB 95**

A polysiloxane based, condensation curing material, exclusively designed for extra oral use in the dental laboratory. Ormalab 95 is useful whenever high hardness and heat-resistance are required.


TECHNICAL DATA:

- Extremely high heat resistance:
- over 120°C.
- Hardness (after 24h): 95(±3) Shore-A.
- Perfect compatibility with cold and heat curing acrylic resins.
- High precision in details reproduction.
- Long lasting dimensional stability.
- High precision in detail reproduction.
- INTENDED USE
- Matrices for fixed and total prosthesis.
- Duplications.
- Prosthesis repair.
- Matrices for temporary crowns and bridges.
- CAD compatibility.

M9500	Ormalab 95 (1600 g / 900 ml)	\$54.38
M9510	Ormalab 95 (5 kg)	\$146.84
M9520	Ormalab 95 (5 kg) + 2 Ormactivator (60 ml)	\$217.49

#### **C-Silicones for Lab**



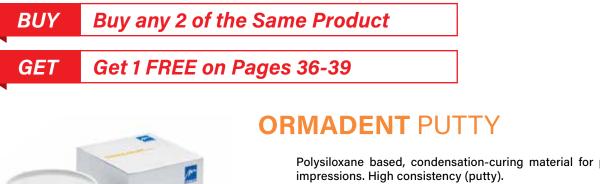

Polysiloxane based, condensation-curing material for preliminary precision

- **TECHNICAL DATA**
- Non-stick material. · Very high compressive strength and recovery from deformation.
- Adequate total working time.
- INTENDED USE Ormadent Putty, due to its excellent initial flowing and final rigidity, guarantees the best results for all impression techniques.

M1000	Ormadent Putty	1500 g (900 ml)	\$64.58



### **ORMACTIVATOR** GEL

Activator Gel for Ormadent Putty and Ormamax Light.

M1090	Ormactivator Gel (35 ml)	\$20.07
M1092	Ormactivator Gel (60 ml)	\$33.97

### **ORMAPLUS** PUTTY (Green and Whte)

Polyvinyl siloxane based, addition curing, elastomeric precision preliminary impression material. Final hardness: 60 Shore-A.

#### INTENDED USE

Base and catalyst pastes mixing ease. Fast intraoral setting time (2'30) which is pleasant for patients. Perfect chemical bonding between putty and light bodied materials. High hydrophilicity and tear resistance. Easy disinfection without dimensional alterations.

Permanent dimensional stability. Perfect compatibility with any kind of plaster (natural or synthetic). Possibility of casting the same impression several times. Odorless and tasteless material for a better patient's comfort.

**TECHNICAL DATA** Fast: Setting Time 4:00 Regular: Setting Time 5:30

M3035	Ormaplus Putty Fast setting: 600 ml / 1060 g (300 ml base paste + 300 ml hardener paste)	\$88.37
M3045	Ormaplus Putty Regular setting: 600 ml / 1060 g (300 ml base paste + 300 ml hardener paste)	\$88.37





### **Condensation Curing Silicone**





## Thermoforming blank made of PETG/EVA, including insulation foil for improved and endured transparency, for the production of:

- occlusion splints for more comfortable wearing characteristics
- hard surface
- soft inner side
- bonds with acrylate clear-transparent

GROUP/ITEM #	DESCRIPTION		SPECIAL
031-D420016	Malibu Thermo Hard/Soft 120 x 2 mm	\$66.41	\$63.25
031-D420016	Malibu Thermo Hard/Soft 120 x 3 mm	\$87.73	\$83.56

#### **BASE PLATE WAX**

Three options available:

**Medium Soft:** dependable wax that boils out clean, reasonably priced and polishes to a high shine.

**Medium Hard:** same as Medium Soft but formulated with higher melting waxes and is specifically designed for warmer climates.

**All Season:** designed to be used most anywhere a good quality, economy wax is needed. This wax boils out clean, polishes nicely, and is formulated to securely hold teeth.

GROUP/ ITEM #	DESCRIPTION	
066-66-013/1LB	Medium Soft size sheets 76mm x 152mm x .060 thick in 1 lb. package.	\$18.20
066-66-014/1LB	Medium Hard size sheets 76mm x 152mm x .060 thick in 1 lb. package.	\$19.34
066-66-012/1LB	All Season size sheets 76mm x 152mm x .060 thick in 1 lb. package.	\$19.76
066-66-016/5LB	Medium Soft size sheets 76mm x 152mm x .060 thick in 5 lb. package.	\$42.47
066-66-017/5LB	Medium Hard size sheets 76mm x 152mm x .060 thick in 5 lb. package.	\$42.47
066-66-015/5LB	All Season size sheets 76mm x 152mm x .060 thick in 5 lb. package.	\$74.18





39

#### AUTOCLAVABLE - IMMEDIATE CUSTOM TRAYS - REUSABLE 24 PIECES

Autoclavable

The Pozzi Super Tray system is designed to produce fine quality and accurate full dentures, and immediate denture impressions. With the Pozzi system, you no longer need to make a custom tray because Pozzi trays are custom trays-ready to use. These trays are easy to handle; you do not need to learn special techniques. You need only apply the procedures that you use with a lab-made custom tray.

Advantages: (1) much more economical than a custom tray (save time and money). (2) Less expensive than disposable trays; these trays are reusable and can be sterilized. (3) Rough retentive surface excludes the use of adhesives. (4) Statistically selected system and packaged in an easy-to-choose kit.

Characteristics: (1) Designed to provide consistent retromolar pad coverage. This prevents errors in applying too much pressure to this area, which is extremely important for the final fitting of the prosthetic appliance. (2) Trays are exceptionally strong and rigid but lightweight and can be modified by trimming for special needs. (3) Can be used with any type of impression material. We suggest the use of silicones like our Major® Ormadent Putty and Ormaplus on page 38 and a very creamy alginate like our Major Alginmax on page 34. (4) Professionally designed with complete anatomic detail.

(5) Hand finished.

Composition: Reusable transparent dark brown polycarbonate.

Technique: No particular technique is suggested. The Pozzi Super Tray System consists of accurate ready-made custom trays.

A caliper is available for measuring the distance between ridges in the patient's mouth or in existing dentures. It helps to expedite proper tray selection.

#### GROUP/ITEM # DESCRIPTION

003-PSTS-K24	Pozzi Super Tray Kit, Kit contains: 15 Lower Trays, 9 Upper Trays & 1 Caliper	\$298.30 / Kit
077-65-0034	Justi Caliper, Replacement	\$19.08



#### Pozzi® Guaranteed Stainless Steel Edentulous Impression Trays

Specially designed trays for complete denture applications. Retention holes and rim are a guarantee for a perfect impression. Available in a kit of 4 upper and 4 lower trays or as single trays.

	GROUP/ITEM # 003-S2806-8	DESCRIPTION Stainless Steel Impression Kit, 8 pcs. \$278.04 Kit contains one each: U61, U62, U63, U64, L66, L67, L68, L69	SPECIAL \$215. <sup>53</sup>
\$34.79ea.	003-2806-U61 003-2806-U62 003-2806-U63 003-2806-U64 003-2806-L66 003-2806-L68 003-2806-L69	Upper Tray, Width 79 mm, Length 70 mm Upper Tray, Width 76 mm, Length 67 mm Upper Tray, Width 70 mm, Length 64 mm Upper Tray, Width 73 mm, Length 60 mm Lower Tray, Width 70 mm, Length 64 mm Lower Tray, Width 73 mm, Length 64 mm Lower Tray, Width 73 mm, Length 67 mm	SPECIAL \$26.91 each



#### **Impression Trays & Accessories**



### Pozzi<sup>®</sup> Guaranteed Standard Stainless Steel Trays



Perforated stainless steel impression trays-the most popular type. Retention is assured by holes over the whole surface of the tray. A rim on the back edge stops the escape of any material, which is most unpleasant for patients. An anatomical design with rounded comfortable rims especially for adaptation to the oral cavity with respect to frenum, jaw tuberosity and retromolar pads. The dimensions come from research in laboratories throughout the world and are able to cover any anatomical situation. This means more accurate and precise impressions.

	GROUP/ITEM # 003-S2801-12	DESCRIPTION Stainless Steel Impression Kit, 12 pcs. \$193.53 Kit contains one each: U1, U2, U3, U4, U5, U6, L1, L2, L3, L4, L5, L6	
\$16.04 ea.	003-2801-U1 003-2801-U2 003-2801-U3 003-2801-U4 003-2801-U5 003-2801-U5 003-2801-L1 003-2801-L2 003-2801-L3 003-2801-L4 003-2801-L5 003-2801-L6	Upper Tray, Width 54 mm, Length 41 mm Upper Tray, Width 61 mm, Length 47 mm Upper Tray, Width 64 mm, Length 50 mm Upper Tray, Width 70 mm, Length 59 mm Upper Tray, Width 75 mm, Length 61 mm Lower Tray, Width 54 mm, Length 41 mm Lower Tray, Width 71 mm, Length 49 mm Lower Tray, Width 71 mm, Length 54 mm Lower Tray, Width 78 mm, Length 54 mm Lower Tray, Width 78 mm, Length 54 mm Lower Tray, Width 78 mm, Length 58 mm	
	003-S2801-6	Stainless Steel Impression Kit, 6 pcs. \$92.16 Kit contains one each: U3, U4, U5, L3, L4, L5 S60.38	
\$16.04 ea.	003-2801-U3 003-2801-U4 003-2801-U5 003-2801-L3 003-2801-L4 003-2801-L5	Upper Tray, Width 64 mm, Length 50 mm Upper Tray, Width 67 mm, Length 59 mm Upper Tray, Width 70 mm, Length 61 mm Lower Tray, Width 71 mm, Length 51 mm Lower Tray, Width 76 mm, Length 54 mm Lower Tray, Width 78 mm, Length 58 mm	

## Pozzi Super-Flexible Mixing Bowls



Their soft flexible texture and sturdy construction make these bowls indispensible for mixing plaster or alginate. They easily conform to the hand, and make mixing a smooth and quick process. These selfcleaning, plasticized polyvinyl bowls are available in four convenient sizes from 350 ml to the max 2250 ml size.

GROUP/ITEM #	DESCRIPTION		
017-5501-10	Flexible Mixing Bowl Small Blue	\$8.74	\$8.49
017-5501-12	Flexible Mixing Bowl Regular Blue	\$9.40	\$9.38
017-5501-14	Flexible Mixing Bowl Large Blue	\$10.58	\$10.45
017-5501-16	Flexible Mixing Bowl Max Blue	\$17.86	\$17.85

### **Massad Low Temp Impression Trays**

#### Massad® Low Temp Denplant Impression Trays - Intro Package

Set of 12 trays, 3 upper and 3 lower sizes

- Clear and strong
- Great for orthodontics, implants, immediate dentures, partial dentures, and crowns
- Fully adjustable in less than 1 minute

#### GROUP/ITEM# DESCRIPTION

difeet /in Elina			SPECIAL
003-689-03	Massad Denplant Tray Intro Kit 12/Pack	\$52.11	\$42.26 each
003-689-01	Massad Edentulous Tray Intro Kit 12/Pack	\$50.59	SPECIAL
003-689-02	Massad Edentulous Tray Combo Kit 12/Pack	\$50.59	\$42.26 each
BUY	Buy any 2 alginates (se	e page	34 and 35)
GET	Get 1 FREE Massad Tra	ay Kit	

DECINI

**Laboratory Specials** 



## **Special Inox Plaster Nippers**

#### Pozzi

#### **THE QUICK DESTROYER**

Stainless steel nippers especially handy for accurate chipping of excess plaster from laboratory work models. Cuts plaster and stone.



2074

GROUP/ITEM # DESCRIPTION 007-444 Special Inox Plaster Nippers \$70.85

SPECIAL \$48.30 ea.

## **Spatulas for Alginate and Plaster**

This line of spatulas was designed by dental engineers for better manipulation of alginate and plaster. The curved tips on items 017-434, 017-435, and 017-436 are the most advanced technical design. They help to mix and become a scoop which is a great help in picking up additional material.



GROUP/ITEM #	DESCRIPTION	
017-430	7 1/4" Spatula	\$17.94
017-431	7″ Spatula	\$17.94
017-432	8 1/4" Spatula	\$18.57
017-433	8" Flexible Spatula	\$18.57
017-434	7 3/4" Flexible Spatula	\$17.94
017-435	7 3/4" Rigid, Heavy Blade Spatula	\$17.94
017-436	8 3/4" Rigid, Heavy Blade Spatula	\$18.57

## Knife & Spatual Kit

A complete kit of 18 different high-quality instruments wrapped in a protective red felt roll



Group/Item# DESCRIPTION 017-4500 Knife and Spatula Kit \$319.53

## GUARANTEE Flask Opening Knife

The dual-purpose flask opening knife is a heavy duty, long-lasting instrument made of carbon steel that extends the length of the handle. The butt end has an extrusion designed to open flasks.



805-487-9868

Offer Expires December 31, 2024-not valid with other offers www.americantooth.com | buy online: store.americantooth.com



## **Plaster Knives**

Autoclavable

### Pozzi Guarantee!

#### High quality and fine craftsmanship set these knives apart.



### **Rigid Crown & Bridge Box**

#### 2" Rigid Crown & Bridge Boxes

Crown and bridge tooth boxes with foam insert. Please call for other available colors. SIZE: 2-1/8" x 1-5/8" x 3/4"



**GROUP/ITEM #** DESCRIPTION 077-R2F-BLU/CL 2" Rigid Crown & Bridge Box w/Foam - Blue 077-R2F-CL/CL 2" Rigid Crown & Bridge Box w/Foam - Clear

\$217.08 \$217.08

CALL TODAY TO PLACE YOUR ORDER

805-487-9868

## Justi<sup>®</sup> Wax Solvent No. 20

An efficient, safe and pleasant-scented solvent used for wax elimination and shining after waxing prosthetic cases. An excellent separator for plaster, stone and any other gypsum product. Great for adding to boil tank water, cleaning counter tops, cleaning tools and equipment, and removing wax from clothing.

070-01-0616

BUY

**GET** 



### Justi<sup>®</sup> Milling & Polishing Cream

An all-purpose cream for milling, finishing and intra-oral polishing of porcelain and acrylic materials.



### **Dental Porcelain Products distributed by ATI**





#### **Paste Stain and Glaze Kit**

The 15-color kit offers A, B, C and D shade stains for customizing shade and 11 basic colors allowing the technician or dentist to create any desired characterization. The kit also includes a 15ml bottle of the Stain and Glaze Medium and a 5g jar of Fluorescent Glaze paste.

GROUP/ITEM #	DESCRIPTION	
111-PSGK16	Paste Stain and Glaze Kit, 16 pcs	\$399.00
1001 1001		
<b>Gum Shade Kit</b> Contains 8 paste one paste liquid.	gum shades, 1 paste	glaze, and

GROUP/ITEM #	DESCRIPTION	
111-PSGKGM	Gum Shade Kit-10 pcs	\$245.00

#### CrystalCeram® Fluorescent Glaze Paste

This is the newest advancement in glazes from the largest U.S. manufacturer of stains and glazes for dental ceramics. Designed specifically for use with Zirconia and Lithium Silicates, CrystalCeram<sup>®</sup> Fluorescent Glaze Paste paints on clear and fires Crystal Clear to help create the most natural esthetic restorations.



**GROUP/ITEM # DESCRIPTION** 111-CCZFG5 CrystalCeram<sup>®</sup> Fluorescent Glaze Paste, 5g **\$50.00** 

#### **Universal Fluorescent Glaze Paste**

CGI Universal Fluorescent Glaze paste is intended for use in glazing all types of dental ceramic restorative materials, especially lithium disilicate, lithium silicate and zirconia restorations. The ideal level of fluorescence enhances the optical properties of the restoration to mimic the appearance of natural teeth under all lighting conditions. The paste permits easy application over the surface of the restoration and saves time during the preparation and application.



GROUP/ITEM #	DESCRIPTION	
111-CFLLFP	Universal Fluorescent Glaze Paste, 5g	\$45.00

#### **Universal Glaze Paste**

CGI Universal Glaze paste is intended for use in glazing all types of dental ceramic restorative materials. Its wide range of firing temperatures allows the technician to adjust to the firing of the specific ceramic to be glazed. The paste permits easy application over the surface of the restoration and saves time during the preparation and application.



GROUP/ITEM #	DESCRIPTION	
111-COGLFP	Universal Glaze Paste, 5g	\$33.00

#### **Universal Glaze Powder**

CGI Universal Low Glaze powder is intended for use in glazing all types of dental ceramic restorative materials. Its wide range of firing temperatures allows the technician to adjust to the firing of the specific ceramic to be glazed.

ú	6				
8					1
1000	-in-	C	G	-	1

GROUP/ITEM #	DESCRIPTION	
111-COGLFS	Universal Glaze Powder, 5g	\$25.00

#### **Universal Fluorescent Glaze Powder**

CGI Fluorescent Low Glaze powder is intended for use in glazing all types of dental ceramic restorative materials, especially lithium disilicate, lithium silicate and zirconia restorations. The ideal level of fluorescence enhances the optical properties of the restoration to mimic the appearance of natural teeth under all lighting conditions.



GROUP/ITEM #	DESCRIPTION	
111-CFLLFS	Universal Fluorescent Glaze Powder, 5g	\$35.00

THE ORIGINAL MASTERPIECE

Bv

*Metal articulator that controls the correct occlusion with maximum precision. No plaster or accessories are needed.* 

#### Featuring Mechanical Fixation of Models Precision Made - Original Type Patented

GALETT®

The Galetti Articulator is solidly constructed of metal alloy casting. It is easy to use, allowing quick fixation of model from chairside to patient. For the dentist, orthodontist and oral surgeon, the Galetti Articulator creates an impressive display that is quick and easy to use for case presentations. It can also be used for study models, as well as for many additional applications for partial and full denture work.

#### **Special Features**

- All anatomical movements possible. Correct occlusion with maximum precision.
- Perfect occlusal-articular balance to prevent cuspidal interference and to obtain a balanced prosthetic appliance.
- Allows easy separation of upper and lower components to form two entirely stable working bases.
- Upper anterior ball-bearing joint allows for free motion of upper appliance to check occlusal-articulator balance.
- Hinge axis with bilateral motion capability for anterior, posterior and lateral movements showing working and balancing movements.
- Adjustable vertical angle pin table to raise or lower the bite.

#### PRECISION INSTRUMENT DO NOT BUY SUBSTITUTES

Group/Item # 001-19893

Galetti Articulator, Complete

\$627.00

SPECIAL

\$501.11 EA



### Articulators

## **Magnett® Articulator**

Features four basic movements: mastication, protrusive, right lateral and left lateral. The protrusive movement can be simulated by placing both thumbs on the upper cross member just inside the condular joints and simply pushing towards the back of the condular springs. To achieve right lateral movement press with the left thumb and to achieve left lateral movement, press with the right thumb against the cross member. Comes with 10 steel mounting plates.

#### **NEW VALUE-PRICED** Addition to the Pozzi Brand

Group/Item # 001-20123 Magnett Articulator **\$134.06** 

Group/Item # 001-20124 Mounting Plates Replacement Parts \$13.12





## **Timeless Articulator**

Justi Timeless Articulator creates an impressive display that is quick and easy to use for case presentations. It can also be used for study models, as well as for many additional applications for partial and full denture work.



#### Free with tooth stock, see page 9 & 11

## **AVM 100M**

A Mid-Range Articulator Without Facebow. A GREAT VALUE!



 Group/Item#
 Description

 080-E3000
 AVM 100M \$210.76

 080-E3020
 Mounting plates 10/pk \$139.36



## Kura Curing Pot with Air Hose



The Justi-Kura<sup>™</sup> Curing Pot has been designed to satisfy the need of the modern technician. It is sturdy, very practical, absolutely safe, and has a large capacity of 6 quarts. The assembling of the pot is done under strict quality control. The price is very competitive considering that it is made entirely of stainless steel. All parts are available individually.

The Justi-Kura<sup>™</sup> Curing Pot is designed for curing Pozzi Ortho Resin, Justi Ortho Resin, and Justi Repair Material.

This system ensures a pourous-free cure. When powder and liquid application is complete, cover appliance with warm water in the pot, move the lid to the closed position and inject air 20 to 30 psi. Allow the resin to cure 15 minutes. Exhaust the air by pressing the pressure release valve until the safety lock drops completely. Open lid, remove appliance and allow to bench cool.

## CHARACTERISTICS & ACCESSORIES

- The Justi-Kura Curing Pot is sturdy, absolutely safe, very practical and has a large capacity of 6 quarts.
- Air Intake Valve Admits compressed air. Valve must be firmly imbedded into the lid.
- Safety Relief Valve Cuts in if the pressure is above the recommended psi and the pressure will be released.
- Pressure Release Valve -To open the pot, press the air release valve until the safety lock drops. Open the lid counterclockwise and lift up.
- Air Pressure Gauge -The pot is preset not to exceed 40 psi.
- Safety Lock Small seal for pressure pot safety lock.
- Air Supply Hose -To inject compressed air into pot.



#### MADE IN THE U.S.A.

GROUP/ITEM #	DESCRIPTION	PRICE
070-46-0001	Justi Kura Pot with air hose	\$470.71
072-46-0002	Air Intake Valve for Justi Kura Pot	\$7.46
072-46-0004	Safety Relief Valve for Justi Kura Pot	\$27.60
072-46-0005	Pressure Release Valve for Justi Kura Pot	\$19.64
072-46-0006	Air Pressure Gauge for Justi Kura Pot	\$44.34
072-46-0007	Small Seal for Pressure Pot Safety Lock	\$5.25
072-28-0023	Air Supply Hose	\$30.20

## **Pozzi Polymerization Flask**

Our Microwave Flask is designed for use during microwave curing procedures. Pozzi's flask for microwave polymerization is made of a reinforced engineered plastic which guarantees excellent physical and mechanical properties.



#### **Benefits**

- Work faster
- Improve the process
- Save money
- Better dimensional stability
- High-temperature resistance
- 100% microwave compatible
- Safe to 2x extertion of pressure
- Plenty of space to work with total upper and lower cases
- Use any Justi microwable resin
- Cleaner than regular procedures

		SPECIAL PRICE	
GROUP/ITEM #	DESCRIPTION		
001-67-0020	Pozzi Polymerization	Flask	\$241.50

## Microflask

Justi is proud to introduce the industrial-strength, two-piece, Super Vipi STG Microflask to be used in the microwave irradiation technique for curing acrylic resins in the fabrication

2 Pieces



- Uniform curing of acrylic resins
- Rapid curing time (13 min)
   Complete curing of acrylic resin with no significant free monomer residual
- Clean removal method without damage to the flask
- Reduced energy costs

of dental prostheses. This revolutionary microwave flask is constructed in 2 sections with 4 stainless steel bolts, sleeves and threads for tightening. Featuring:

#### GROUP/ITEM # DESCRIPTION

(	071-67-0004	
C	071-67-0006	
(	080-67-0007	
(	071-67-0008	

#### Justi Super Vipi STG Microflask \$228.85 STG System Instructional Video \$46.67 2 Replacement Bolts \$7.92 Allen Wrench for VP STG System \$3.97

#### **IMPORTANT NOTICE**

Pack denture with the standard method trial pack using a manual or pneumatic press to a maximum of 200 PSI. Do not over tighten bolts. Tap flask with a wooden mallet only.

Offer Expires December 31, 2024–not valid with other offers www.americantooth.com | buy online: store.americantooth.com





## Pozzi<sup>®</sup> Acrylic Trimming Bur & Laboratory Carbide Burs

Acrylic Trimming Bur: Specially designed for rapid reduction on acrylics. Cuts faster and stays sharper longer than any other conventional bur on the market. Will not clog. Taper-fissure type-one

size only.

Laboratory Carbide Burs: For trimming, carving and finishing acrylic appliances. These carbide burs provide the ultimate in toughness and durability. Precision-ground grooves are razor sharp and assure resistance to gumming and clogging. All carbide burs, except Inverted Cone, feature crosscut, chip-breaker blades. These carbide burs are now coated with gold titanium nitride (GTN) which makes them more durable.



## **Pozzi® Carbide Burs for Bench Motors**



Coated with gold titanium nitride for durability.

Unique designs of these bench burs allow for faster finishing of acrylics and leaves a very smooth surface. Top-quality, long-lasting, fast-cutting, dual-cut carbide burs; finely tooled for maximum efficiency.

The most popular acrylic finishing burs. Cuts so clean you reduce both pumicing and polishing to a minimum.

GROUP/ITEM #	DESCRIPTION	
013-19796/M	6 mm Diam. Carbide Bur 3/8″ Width Head Egg	\$78.15
013-19796/L	6 mm Diam. Carbide Bur 1/2″ Width Head Egg	\$102.72
013-19797/M	6 mm Diam. Carbide Bur 3/8" Width Head Taper	\$77.12
013-715G88AD	88A Carbide Bur GT-COAT 3/32"	\$30.03




## Pages 50-56

### Buy 2 Boxes Same size Get 1 FREE SAME SIZE

### Malibu-Pozzi Products - Non Contaminating - Compressed Abrasive Formula



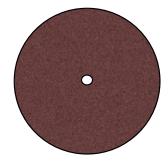


008-39202 PC1

008-39202 PC2

#### Pin Cutters for Chrome-Cobalt & Non-Precious Alloys

Series PC - Available in 35mm diameter, 0.7mm thickness; and 40mm x 1mm thickness.

GROU	JP/ITEM #	DESCRIPTION	
008-3	9202 PC1	Malibu Pin Cutters 35mm x 0.7mm 100/pk.	\$26.29
008-3	9202 PC2	Malibu Pin Cutters 40mm x 1mm 100/pk.	\$33.91

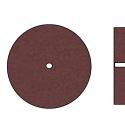

008-39102 AW

#### Abrasive Wheels for Chrome-Cobalt & Non-Precious Alloys

Series AW - available in 22mm diameter, 3.5mm thickness.

GROUP/ITEM #	DESCRIPTION	
008-39102 AW	Malibu Wheel 22mm x 3.5mm 100/pk.	\$30.40



008-39302 RD1



008-39302 RD2



008-39402 S0



008-39402 S1



008-39402 S2

#### **Roughing Discs for Chrome-Cobalt & Non-Precious Alloys**

Series RD - available in 35mm diameter x 1.7mm thickness; and 32mm diameter x 2.2mm Thickness.

GROUP/ITEM #	DESCRIPTION	
008-39302 RD1	Malibu Roughing Discs 35mm x 1.7mm 100/pk.	\$38.41
008-39302 RD2	Malibu Roughing Discs 32mm x 2.2mm 100/pk.	\$38.41

#### Separators for Chrome-Cobalt & Non-Precious Alloys


Series S - Available in 22mm diameter, 0.2mm and 0.6mm thickness; 25mm diameter and 0.7mm thickness.

GROUP/ITEM #	DESCRIPTION	
008-39402 S0	Malibu Separators 22mm x 0.2mm 100/pk.	\$46.69
008-39402 S1	Malibu Separators 22mm x 0.6mm 100/pk.	\$17.09
008-39402 S2	Malibu Separators 25mm x 0.7mm 100/pk.	\$20.47



### **Compressed Aluminum Oxide**

#### Mounted Stones for Chrome-Cobalt Alloy & Non-Precious Alloys

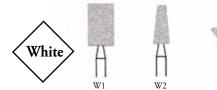



#### White Color with HP type shank.

For gold and porcelain. Made with compressed aluminum oxide.

GROUP/ITEM #	DESCRIPTION	
008-36085 W1, W2, W3	Series W Mounted Stones Pkg. 72/pkg.*	\$26.88



#### **Mounted Stones for Precious Alloys**

Series G mounted stones. Light green color. HP type shank. For gold. Compressed aluminum oxide. Non-contaminating.

GROUP/ITEM #	DESCRIPTION	
008-35048	Mounted Stones for Precious Alloys Pkg. 12 pcs.*	\$15.84



#### **Mounted Stones for Porcelain**

Series P mounted stones. Dark Green color. HP type shank. Non-contaminating.

GROUP/ITEM #	DESCRIPTION	
008-36084	Mounted Stones for Porcelain Pkg. 12 pcs.*	\$23.05

V1	V2	V3	V6
Mounted Stones f Series V mounted stone		reen. HP type	shank. Non-contaminating.


GROUP/ITEM #	DESCRIPTION	
008-34084	Mounted Stones for Resin. Pkg. 6 pcs.*	\$16.62

\*Please specify size when ordering on all W, G, V, & P stones.

### **Mounted Stones for Chrome-Cobalt Alloy & Non-Precious Alloys**


Series M mounted stones for chrome-cobalt alloy. Available in the shapes pictured. Blue in color. USA specification. Medium type, cuts fast. Shank 2.35 mm. diameter. Made of compressed aluminum oxide. Non-contaminating.

Guaranteed to cut fast. Recommended speed is between 20,000 and 28,000 rpm.

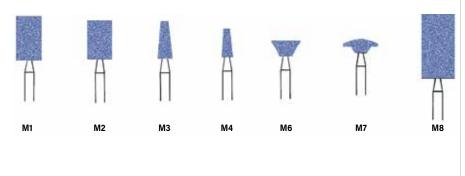

GROUP/ITEM #	DESCRIPTION	
008-32084	Mounted stones for Chrome-cobalt alloy. Pkg. 100 pcs. Series M Please specify size when ordering on all M stones.	\$37.43

## **Abrasive Wheels for Porcelain & Gold**

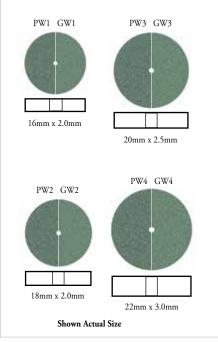

Series PW - Compressed aluminum oxide abrasive wheels for	
porcelain. Non-contaminating.	

Series GW - Compressed aluminum oxide abrasive wheels for precious alloys and gold. Non-contaminating.

GROUP/ITEM #	DESCRIPTION	
008-39084 PW1	Abrasive Wheels for Porcelain, 16mm x 2.0mm, Pkg. 12 pcs.	\$14.43
008-39084 PW2	Abrasive Wheels for Porcelain, 18mm x 2.0mm, Pkg. 12 pcs.	\$14.43
008-39084 PW3	Wheels for Porcelain, 20mm x 2.5mm, Pkg. 12 pcs.	\$14.43
008-39084 PW4	Abrasive Wheels for Porcelain, 22mm x 3.0mm, Pkg. 12 pcs.	\$14.43
008-38584 GW1	Abrasive Wheels for Precious Alloys and Gold, 16mm x 2.0mm, Pkg. 12 pcs.	\$14.43
008-38584 GW2	Abrasive Wheels for Precious Alloys and Gold, 18mm x 2.0mm, Pkg. 12 pcs.	\$14.43
008-38584 GW3	Abrasive Wheels for Precious Alloys and Gold, 20mm x 2.5mm, Pkg. 12 pcs.	\$14.43
008-38584 GW4	Abrasive Wheels for Precious Alloys and Gold, 22mm x 3.0mm, Pkg. 12 pcs.	\$14.43

### Pozzi Abrasives - "Mario" Series for use on Chrome Cobalt Alloy - Made in the U.S.A.


Thick



**Cut-off Wheel - Large / Thick** (11/2" x .040") (38.1mm x 1 mm

GROUP/ITEM #	DESCRIPTION	
010-19921/210	Cut-off Wheel - Large / Thick Pkg. 100 pcs.	\$27.15



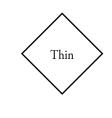

Cut-off Wheel - Large / Thin (11/2" x .025") (38.1mm x 0.63 mm)

GROUP/ITEM #	DESCRIPTION	
010-19922/211	Cut-off Wheel - Large / Thin Pkg. 100 pcs.	\$27.15

#### Pozzi Abrasives - "Mario" Series for use on Chrome Cobalt Alloy - Made in the U.S.A.



#### Cut-off Wheel - Small/Thick

(11/4" x 1/16") (31.8mm x 1.6mm)



High Speed Separating Disc New Improved Red

(15/16" x .025") (23.8mm x 0.63mm)



High Speed Separating Disc - Regular Brown

(15/16" x .025") (23.8mm x 0.63mm)

GROUP/ITEM #	DESCRIPTION		GROUP/ITEM #	DESCRIPTION		GROUP/ITEM #	DESCRIPTION	
010-19923/212	Cut-off Wheel - Small/ Thick Pkg. 100 pcs.	\$27.15	010-19924/221	High Speed Separating Disc Pkg. 100 pcs.	\$11.55	010-19925/222	High Speed Separating Disc Pkg. 100 pcs.	\$10.41



DESCRIPTION

Disc Pkg. 100 pcs.

High Speed Separating

**Regular Double Cutting - Grey** 

(7/8" x .023") (22.2mm x 0.6mm)

Separating Disc

**GROUP/ITEM #** 

010-19926/223



\$12.52





Separating Disc - Regular Safe Side - Grey (7/8" x.023") (22.2mm x 0.6mm)

//0	x .023	) (22.2000 x 0.00000)	

GROUP/ITEM #	DESCRIPTION	
010-19927/224	Separating Disc Pkg. 100 pcs.	\$12.52

### Pozzi Silicone Abrasives - For Amalgams, Silicates, Composites & Precious Alloys

Multi Use         Light Grey         Special silicone abrasives for grinding - light grey - Pkg. 12 pcs.         GROUP/ITEM #       026-19994         Multi Use	PUIW	PU2H PU2W	PU3H PU3W	PU4W	PU5H PU5W
Violet Special silicone abrasives for polishing - violet - Pkg. 12 pcs. GROUP/ITEM # 026-19995 \$31.49 *WHEN ORDERING SPECIFY: For example PK1W means only available in contra-angle. H - For Straight Handpieces W - For Contra-angle Handpieces	Â	PG2H	РСЗН		PG5H
0 2.35	PG1W	PG2W	PG3W	PG4W	PG5W



#### **Pozzi® Abrasive Universal Series**

### Pozzi Rubber Abrasives "Mario" Series - For use on Chrome Cobalt Alloy

O Rubber Point -	Grey	Pul		pecial Red		Blue Magic Rubber Point -	Rhue Magic	
(1/4" x 7/8") (6.4mm x 22.2mm)			•		(1/4" x 7/8") (6.4mm x 22.2mm)			
GROUP/ITEM #	DESCRIPTION	GR	OUP/ITEM #	DESCRIPTION		GROUP/ITEM #	DESCRIPTION	
010-19915/203	Rubber Point - Grey Pkg. 100 pcs.	\$16.36 010	)-19917/205	Rubber Point - Special red Pkg. 100 pcs.	\$17.15	010-19919/207	Rubber Point - Blue Magic Pkg. 100 pcs.	\$17.15
Rubber Wheel • 5/8") (15.9mm)	- Knife Edge							
GROUP/ITEM #	DESCRIPTION			• Knife	$\mathbf{>}$			
010-19920/208	Rubber Wheel - Knife E Pkg. 100 pcs.	dge \$26.59		Edge			Made in U.S.A.	

### Pozzi Rubber Abrasives - "Mario"



### Pozzi Abrasives - "Mario" Series - For use on Porcelain & Precious & Semi-Precious Metals

Ceramic Alloy Finishing Wheel (7/8" x 1/8") (22.2mm x 3.2mm)	Porcelain Grinding Wheel (7/8" x 1/8") (22.2mm x 3.2mm)	<b>Porcelain Grinding Wheel</b> (3/4" x 3/32") (22.2mm x 3.2mm)
GROUP/ITEM # 011-19928/301 \$13.87	GROUP/ITEM # 011-19929/3021 \$13.87	GROUP/ITEM # 010-19930/3022 \$13.87



### Pozzi Abrasives - "Mario" Series - For use on Porcelain & Precious & Semi-Precious Metals



**Porcelain Detail Grinding Wheel** 



Fine Porcelain Polishing Wheel




#### **Fine Porcelain Polishing Wheel**


(5/8" x 1/16") (15.9mm x 1.6mm) – Pkg. 25 pcs.	(7/8" x 1/8") (22.2mm x 3.2mm) – Pkg. 25 pcs	(5/8" x 1/16") (15.9mm x 1.6mm) - Pkg. 25 pcs.			
GROUP/ITEM # 011-19931/3023 \$13.87	GROUP/ITEM # 011-19932/3031 \$13.87	GROUP/ITEM # 011-19933/3032 \$12.71			
•					
Extra Fine Porcelain Polishing Wheel	Extra Fine Porcelain Polishing Wheel	Knife-Edge Extra Fine Porcelain Polishing Wheel			
7/8" x 1/8") (22.2mm x 3.2mm) – Pkg. 25 pcs.	(5/8" x 1/16") (15.9mm x 1.6mm) - Pkg. 25 pcs.	(5/8″) (15.9mm) – Pkg. 25 pcs			
GROUP/ITEM # 011-19934/3041 \$12.71	GROUP/ITEM # 011-19935/3042 \$12.71	<b>GROUP/ITEM #</b> 011-19936/304K \$12.71			
	A 1				
( o ) 🗆		$(\circ)$			
	Cup Shape Separating Disc,				
White Gloss Polishing Wheel	Double Cutting (7/8" x .023") – Pkg. 25 pcs.	Extra Thin Separating Disc			
7/8" x 1/8") (22.2mm x 3.2mm) – Pkg. 25 pcs.		(7/8" x .015") (22.2mm x 0.38mm) – Pkg. 25 pcs.			
GROUP/ITEM # 011-19937/305 \$5.96	<b>GROUP/ITEM #</b> 011-19938/306 \$5.96	<b>GROUP/ITEM #</b> 011-19939/307 \$5.96			
•	•				
Razor Thin Separating Disc	Slim Disc Separating Di	sc			
(7/8" x .009") (22.2mm x 0.22mm) Pkg. 25 pcs					
<b>GROUP/ITEM #</b> 011-19940/308	\$17.15 GROUP/ITEM # 011-19941/	309 \$17.15			



### Pozzi Silicone Carbide Grinding & - Polishing Wheels for Multiple Uses

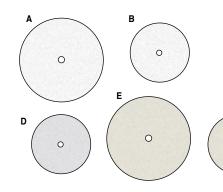

- A Porcelain Grinding Wheel, 7/8" x 1/8"
- B Porcelain Detail Grinding Wheel, 5/8" x 1/16"
- C Fine Porcelain Polishing Wheel, 7/8" x 1/8"

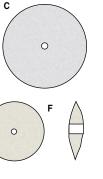

- D Fine Porcelain Polishing Wheel, 5/8" X 1/16"
- E Extra Fine Porcelain Polishing Wheel, 7/8" x 1/8"
- F Extra Fine Porcelain Polishing Wheel, 5/8" Knife Edge

### Multi Use

USE: For polishing natural and artificial teeth, filling materials, acrylics, porcelain, stainless steel and precious metals. Suitable for use up to maximum of 18,000 rpm. When working at low speeds the disc is quite flexible and can be used to polish and finish any inaccessible areas. When working at higher speeds, the discs naturally become less flexible and are used on larger occlusal and labial surfaces. Dry and wet polishing possible.

GROUP/ITEM #	DESCRIPTION	
009-19839*	Porcelain Grinding Wheels, pkg. 12 pcs, Type A, B, C, D, E or F. One type per pack- age, Specify type when ordering.	\$11.56


\*example 19839 A


### **Pozzi Silicone Rubber Abrasives**


#### **Pozzi Silicone Rubber Abrasives**



#### Pozzi Silicone Abrasives - For Ceramic & Porcelain

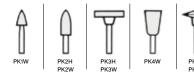

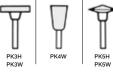

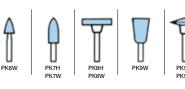





#### Light Grey

GROUP/ITEM

Cepoline SG special silicone abrasives for grinding - light grey. - Pkg. 12 pcs.

B	I	u	e

Cepoline F special silicone abrasives for detail grinding blue. - Pkg. 12 pcs.

#### Grev

Cepoline Lux special silicone abrasives for polishing grey - Pkg. 12 pcs.

EM #	026-19993	\$31.49	GROUP/ITEM #	026-19993	\$31.49	GROUP/ITEM #	026-19993	\$31.49

Offer Expires December 31, 2024-not valid with other offers www.americantooth.com | buy online: store.americantooth.com



The Pozzi stainless steel crown line is manufactured to meet the needs of all facets of dentistry. The posterior molar crowns are designed with a shallow occlusal anatomy, allowing for more healthy tooth structure after tooth reduction. The shallow occlusal cusps permit complete compatibility with opposing teeth and reduced cuspal interference during normal mastication. The crowns are parallel on the mesial and distal walls. Unlike bell-shaped crowns, these crowns will not rock or cause separation when two or more adjoining crowns are used. We have made these crowns with a thickness of 0.18 mm that prevents occlusal wear-through that sometimes occurs with the use of thinner crowns. Unlike some of the other crowns made with thinner and less durable

materials, our crowns may be welded too. These are the crowns you want for making space maintainers or Herbst appliances, or for a semipermanent restoration. Through our own research, we have determined these crowns to be of the ideal length: short enough to eliminate excess crown trimming, and long enough to provide optimum length and contour. To make our line complete to the user, we offer first and second primary molars, first and second permanent molars and first and second bicuspids. Pozzi Style 2000 crowns are a match for Unitek<sup>®</sup> Crowns. The Pozzi<sup>®</sup> Crown numbers and sizes are the same as Unitek<sup>®</sup>. Unitek<sup>®</sup> is not a trademark of American Tooth Industries.







#### 1<sup>st</sup> Primary Molar UR 1 GROUP/ITEM # DESCRIPTION 2 3 First Primary Molar Crown Kit 040-40-0022 4 72 pcs., 7 sizes 5 Reorder pack of 5. 6



#### 4 5 6 7 7 LR LL 1 1 2 2 З 3 4 4 5 5 6 6

7

UR

1

UL

1

2

3

7

UL

1

2<sup>nd</sup> Primary Molar

GROUP/ITEM #	DESCRIPTION
040-40-0023	Second Primary Molar Crown Kit 72 pcs., 7 sizes
040-40-23 (+size)	Reorder pack of 5. Specify crown size and quadrant. *Example: 40-23UL1
MOLAR CROWN KIT SPECIAL: \$283.76	CROWN SPECIAL: BUY 5 PKS, PAY \$16.85 PER PK, MIX & MATCH
Regular kit price \$302.04	Reorder a pack of five crowns \$20,60

1<sup>st</sup> Permanent Molar

GROUP

040-40-

040-40-2

/ITEM #	DESCRIPTION
0024	First Permanent Molar Crown Kit 72 pcs., 8 sizes
24 (+size)	Reorder pack of 5. Specify crown size and quadrant. *Example: 40-24UL-1



Reorder a pack of five crowns \$20.60 UR UL 1 1 2 2 3 3

4

5

6

7

1

2

З

4

5

6

7

8

LL

4

5

6

7

LR 1

2

3

4

5

6

7

8

КΠ

#### Pozzi<sup>®</sup> Stainless Steel Crowns







### 2<sup>nd</sup> Permanent Molar

		2
GROUP/ITEM #	DESCRIPTI	
040-40-0025	Second Permanent Molar Crown Kit - 72 pe	cs., 4
040-40-0025	8 sizes	5
	Reorder pack of 5.	6
040-40-25(+size)*		· · · ·
	*Example: 40-25UL-1	LR
		1
MOLAR CROV	CROWN SPECIAL:	2
KIT	BUY 5 PKS.	3
SPECIAL:	PAY \$16.85 PER PK,	5
\$283.76	MIX & MATCH	6
Regular kit price \$30	2.04 Reorder a pack of	7
	five crowns \$20.60	8
St Bicus		<u>UR</u> 00 0
ROUP/ITEM #	·	00
ROUP/ITEM #	DESCRIPTION First Permanent Bicuspid Crown	00 0 1 2 3
ROUP/ITEM #	DESCRIPTION First Permanent Bicuspid Crown Kit - 72 pcs., 8 sizes (uppers) 7 sizes (lowers) Reorder pack of 5.	00 0 1 2 3 4
ROUP/ITEM # 0-40-0020	DESCRIPTION First Permanent Bicuspid Crown Kit - 72 pcs., 8 sizes (uppers) 7 sizes (lowers) Reorder pack of 5. Specify crown size and	00 0 1 2 3 4 5
ROUP/ITEM #	DESCRIPTION First Permanent Bicuspid Crown Kit - 72 pcs., 8 sizes (uppers) 7 sizes (lowers) Reorder pack of 5. Specify crown size and quadrant.	00 0 1 2 3 4 5 6
ROUP/ITEM #	DESCRIPTION First Permanent Bicuspid Crown Kit - 72 pcs., 8 sizes (uppers) 7 sizes (lowers) Reorder pack of 5. Specify crown size and	00 0 1 2 3 4 5 6 <b>LR</b>
ROUP/ITEM #	DESCRIPTION First Permanent Bicuspid Crown Kit - 72 pcs., 8 sizes (uppers) 7 sizes (lowers) Reorder pack of 5. Specify crown size and quadrant. *Example: 40-20UL1	00 0 1 2 3 4 5 6 <b>LR</b> 0
ROUP/ITEM # 10-40-0020 10-40-20 (+size)* MOLAR CROV	DESCRIPTION First Permanent Bicuspid Crown Kit - 72 pcs., 8 sizes (uppers) 7 sizes (lowers) Reorder pack of 5. Specify crown size and quadrant. *Example: 40-20UL1 WN CROWN SPECIAL: BUY 5 PKS,	00 0 1 2 3 4 5 6 <b>LR</b> 0 1
ROUP/ITEM # 10-40-0020 10-40-20 (+size)* MOLAR CROV KIT SPECIAL	DESCRIPTION First Permanent Bicuspid Crown Kit - 72 pcs., 8 sizes (uppers) 7 sizes (lowers) Reorder pack of 5. Specify crown size and quadrant. *Example: 40-20UL1 WN CROWN SPECIAL: BUY 5 PKS, PAY \$16.85 PER PK,	00 0 1 2 3 4 5 6 <b>LR</b> 0
ROUP/ITEM # 10-40-0020 10-40-20 (+size)* MOLAR CROV KIT SPECIAL: \$283.76	DESCRIPTION First Permanent Bicuspid Crown Kit - 72 pcs., 8 sizes (uppers) 7 sizes (lowers) Reorder pack of 5. Specify crown size and quadrant. *Example: 40-20UL1 WN CROWN SPECIAL: BUY 5 PKS, PAY \$16.85 PER PK, MIX & MATCH	00 0 1 2 3 4 5 6 <b>LR</b> 0 1 2
ROUP/ITEM # 10-40-0020 10-40-20 (+size)* MOLAR CROV KIT SPECIAL: \$283.76	DESCRIPTION         First Permanent Bicuspid Crown         Kit - 72 pcs., 8 sizes (uppers)         7 sizes (lowers)         Reorder pack of 5.         Specify crown size and         quadrant.         *Example: 40-20UL1         VN         CROWN SPECIAL:         BUY 5 PKS,         PAY \$16.85 PER PK,         MIX & MATCH         2.04	00 0 1 2 3 4 5 6 <b>LR</b> 0 1 2 3
ROUP/ITEM # 40-40-0020 40-40-20 (+size)* MOLAR CROV KIT SPECIAL: \$283.76	DESCRIPTION First Permanent Bicuspid Crown Kit - 72 pcs., 8 sizes (uppers) 7 sizes (lowers) Reorder pack of 5. Specify crown size and quadrant. *Example: 40-20UL1 WN CROWN SPECIAL: BUY 5 PKS, PAY \$16.85 PER PK, MIX & MATCH	00 0 1 2 3 4 5 6 <b>LR</b> 0 1 2 3 4
ROUP/ITEM # 40-40-0020 40-40-20 (+size)* MOLAR CROV KIT SPECIAL:	DESCRIPTION         First Permanent Bicuspid Crown         Kit - 72 pcs., 8 sizes (uppers)         7 sizes (lowers)         Reorder pack of 5.         Specify crown size and         quadrant.         *Example: 40-20UL1         VN         CROWN SPECIAL:         BUY 5 PKS,         PAY \$16.85 PER PK,         MIX & MATCH         12.04	00 0 1 2 3 4 5 6 <b>LR</b> 0 1 2 3 4 5

UR

lui

	•		UL
GROUP/ITEM #	DESCRIPTION	0	0
040-40-0021	Second Permanent Bicuspid Crown Kit - 72 pcs., 8 sizes (uppers) 7 sizes (lowers)	1 2 3 4 5	1 2 3 4 5
040-40-21(+size)*	Reorder pack of 5. Specify crown size and quadrant. *Example: 40-21UL1	6 7 <b>LR</b>	6 7 LL
MOLAR CROV KIT SPECIAL: \$283.76	VN CROWN SPECIAL: BUY 5 PKS, PAY \$16.85 PER PK, MIX & MATCH	0 1 2 3 4	0 1 2 3 4
Regular kit price \$30		5	5

### Pozzi<sup>®</sup> **Space Maintainers**

These preformed components are available to the laboratory that chooses to form a custom band or crown space maintainer. You simply solder the tube to the band/crown, insert the maintainer rest into the tubing, and crimp it at the desired point. Your doctor can even do the crimping in the office, thereby eliminating any chance of misfitting.



### Pozzi® straight length and spool wires

#### Pozzi<sup>®</sup> Retainer Wire (Medical Specification)

Diamond drawn: yields a perfectly round wire.

- Medium hard temper: a true 1/2 hard wire.
- Easily formed: ideal for hand forming.
- Non-brittle: can be bent 12-14 times before breaking to prevent remakes.





#### DESCRIPTION

015-19907-2

015-19907-3

015-19907-4

candard Wire

 Retainer Wire Straight lengths of 14", tube .025 10 pcs.
 \$4.92

 Retainer Wire Straight lengths of 14", tube .028 10 pcs.
 \$4.92

 Retainer Wire Straight lengths of 14", tube .030 10 pcs.
 \$5.41

 Retainer Wire Straight lengths of 14", tube .030 10 pcs.
 \$5.45

 Retainer Wire Straight lengths of 14", tube .032 10 pcs.
 \$5.45

 Retainer Wire Straight lengths of 14", tube .036 10 pcs.
 \$5.45

BUY 4 SIZES AND GET 1 SIZE FREE OF EQUAL OR LESSER VALUE

Standard Round Wire Size: .032 10/pkg \$5.61

Standard Round Wire Size: .036 10/pkg \$5.61

Standard Round Wire Size: .040 10/pkg \$5.61 Standard Round Wire Size: .045 10/pkg \$5.61

Standard Round Wire Size: .051 10/pkg \$5.61

### Pozzi® Standard Round Wire

(Medical Specification)

- High Elasticity, which permits greater manipulation of this 3/4 hard wire without danger of cracking or breaking.
- High yield point which gives the wire a greater resistance to permanent deformation.
- Uniform mechanical properties from wire to wire and batch to batch.
- Packaged in 10 14" straight lengths per tube. When ordering please specify item number and size.

Standard Bound Wire Size: .019 10/pkg. \$4.70

GROUP/ITEM #	DESCRIPTION
015-19908-05	Standard Round Wire Size: .012 10/pkg. \$4.70
015-19908-07	Standard Round Wire Size: .014 10/pkg. \$4.70
015-19908-08	Standard Round Wire Size: .015 10/pkg. \$4.70
015-19908-09	Standard Round Wire Size: .016 10/pkg. \$4.70
015-19908-10	Standard Round Wire Size: .017 10/pkg. \$4.70
015-19908-11	Standard Round Wire Size: .018 10/pkg. \$4.70

#### Pozzi<sup>®</sup> Cobalt Clasp Crozat Wire (Medical Specification)

- Chromium cobalt alloy. Specifically designed for the Crozat method and for fabrication of the OSAYS appliance.
- · Heat-treatable: May be soldered in open flame.
- Available in 15-23 gauge.

015-19908-12

 Coiled spools. A softer, more malleable wire saves time when making intricate bends.



015-19908-13

015-19908-14

015-19908-15

015-19908-16

015-19908-17

015-19908-18

015-19908-19

015-19908-20

#### Pozzi<sup>®</sup> Standard Round Wire (Medical Specification)

015-19908-21

015-19908-22

015-19908-23

015-20016-1

015-20016-2

- 3/4 hard wire for those who want more spring to their wire.
- Surgical-grade stainless steel.
  Spooled wire allows for increased
- tensile strength and less waste.



#### SPECIAL: BUY 3 PIECES SAME SIZE AND GET 4<sup>th</sup> PIECE FREE (ANY SIZE - SAME PRODUCT ONLY, OF EQUAL OR LESSER VALUE)

Standard Round Wire Size: .020 10/pkg. \$4.80

Standard Round Wire Size: .021 10/pkg. \$4.80

Standard Round Wire Size: .022 10/pkg. \$4.80

Standard Round Wire Size: .024 10/pkg. \$4.80

Standard Round Wire Size: .025 10/pkg. \$4.80

Standard Round Wire Size: .026 10/pkg. \$5.61

Standard Round Wire Size: .028 10/pkg. \$5.61

Standard Round Wire Size: .030 10/pkg. \$5.61

Gauge Size	mm.	inch	ozs.#	ft.	Group/Item#	
22/23	0.6	.024	2/3	22	015-410/6	\$51.31
21	0.7	.028	3/4	22	015-410/7	\$51.31
20	0.8	.032	1	22	015-410/8	\$51.31
19	0.9	.036	1 1/3	22	015-410/9	\$51.31
18	1.0	.040	3/4	11	015-410/1.0	\$56.44
17	1.1	.044	1	11	015-410/1.1	\$56.44
16	1.2	.048	1 1/4	11	015-410/1.2	\$56.44
15	1.3	.052	1 1/3	11	015-410/1.3	\$56.44

BEST IN THE WORLD	



Gauge						
Size	mm.	inch	OZS.	# ft.	Group/Item	
# 21	0.7	.028	3	90	015-400/7	\$32.60
20	0.8	.032	3	68	015-400/8	\$32.60
19	0.9	.036	3	54	015-400/9	\$32.60 \$32.60
18	1.0	.040	3	44	015-400/1.0	\$32.60
17	1.1	.044	3	36	015-400/1.1	\$32.60
16	1.2	.048	3	30	015-400/1.2	\$32.60
15	1.3	.052	3	26	015-400/1.3	\$32.60



59

#### Pozzi Super Ball Clasps - Series 2000

This improved ball clasp is made of high-grade #304 stainless steel wire and accepts severe bends without fracturing. The round head is approximately twice the diameter of the wire shaft and will not irritate tissue. These clasps are made from diamond-drawn wire and have a bright, clean finish. In our opinion, other ball clasps on the market today have a lower tensile strength and an oval-shaped head that does not possess adequate retention required by retainers. If you are looking for a better ball clasp to give superior retention on your removable appliances, try the 2000 series.

Available in three sizes.						
GROUP/ITEM #	DESCRIPTION					
013-19814	Ball Clasp 2000, Size .028, Pkg. 100 pcs.	\$38.78				
013-19815	Ball Clasp 2000, Size .032, Pkg. 100 pcs.	\$38.78				
013-19815/6	Ball Clasp 2000, Size .036, Pkg. 100 pcs.	\$38.78				







#### **Pozziloy** (Medical Specification)

Pozziloy retainer wire, is our own special formulation of nickel cobalt wire, available in a blue, soft, heat-treatable (1400/1600 N/sg.mm) tensile. This chromium-cobalt base alloy is highly stainless and mechanically more resistant to breakage. Each wire is color coded blue at the tip, and is 14 inches in length. This wire is made from a vacuum-fused alloy, diamond-drawn to obtain leveled and highly resistant forces. Packaged in a ten pack, it is available in .028, .030, .032 and .036.

N	riloy Retainer Wire Addical Specification a Non- 918 Physics Pre-		,
GROUP/ITEM 040-40-0002 040-40-0003 040-40-0004 040-40-0005	<ul> <li># DESCRIPTION</li> <li>Pozziloy .028 –</li> <li>Pozziloy .030</li> <li>Pozziloy .032 _</li> <li>Pozziloy .036</li> </ul>	\$28.51 ea. \$34.18 ea.	SPECIAL: BUY 4 GET 1 FREE! OF EQUAL OR LESSER VALUE MIX & MATCH

### Pozzi<sup>®</sup> Wire & Sheet Metal Gauge

#### **Precision Instrument**

Light-weight stainless steel. Measures round and rectangular wire with gauge and inch sizes. The only device which helps to verify sizes and maintains guality control. The wire gauge is a template which is used to measure the gauge or thickness of sheet metal or wire. Around the edge are graduated, numbered slots which fit wire or sheet metal. The metal or wire is tried in the various slots until it just fits without being forced.

GROUP/ITEM # DESCRIPTION Wire Gauge \$60.78 014-19895



### **Soldering Tweezers**

Stainless steel soldering tweezers with plastic handles for protection from heat of soldering. Available in two styles: straight or curved tipped.

#### GROUP/ITEM # DESCRIPTION

007-6400-1	Soldering Tweezers Straight Tip 6½	\$17.17
007-6400-2	Soldering Tweezers Curved Tip 61/4	\$17.17



Pozzi All Purpose	e Scissors
7" long, stainless steel, autoclavable 143 <sup>0</sup> C, with plastic handles.	SPECIAL 620.05

GROUP/ITEM # DESCRIPTION 007-20012 All Purpose Scissors \$37.68





Half round, semi-hard, type-304 tempered stainless steel wire to be used for reinforcement of dentures and prosthodontic appliances. Also used to fabricate special clasps. Precision cut to 12" (30.5 cm) in length. Diamond drawn for brightness, ultrasonically cleaned for dental use. Available in six (6) different sizes. Packaged in 6 pieces of 12" straight length per tube.

GROUP/ITEM # DESCRIPTION

015-20074	12" Half Round Wire .90 mm 6/pk	\$5.17 ea.	
015-20075	12" Half Round Wire 1 mm 6/pk	\$5.17 ea.	
015-20076	12" Half Round Wire 1.2 mm 6/pk	\$5.17 ea.	
015-20077	12" Half Round Wire 1.25 mm 6/pk	\$7.47 ea.	ľ
015-20078	12" Half Round Wire 1.5 mm 6/pk	\$7.87 ea.	
015-20079	12" Half Round Wire 1.7 mm 6/pk	\$8.91 ea.	

#### **JAL: BUY** 21 QUAL OR LESSER VALUE MIX & MATCH


### **Pozzi® Lingual Bars**

Stainless steel lingual bars in four sizes: small, medium, large, extra large. More retention due to large serrated anchorage surface.

GROUP/ITEM # DESCRIPTION

015-20081 Lingual Bar, Small - 55 mm. 1/pkg. Lingual Bar, Medium - 65 mm. 1/pkg. 015-20082 \$3.04 ea. 015-20083 Lingual Bar, Large - 75 mm. 1/pkg. Lingual Bar, X-Large - 85 mm. 1/pkg. 015-20084

MIX&MATC

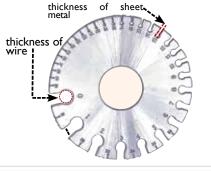



Pozzi<sup>®</sup> Wire/Allov





#### **Premium Ticonium Wire**

Premium quality, round clasp wire specially suited for casting with PREMIUM 100 alloy. 12" length (30cm) each wire. Comes 5 Sticks in a tube.

GROUP/ITEM #	DESCRIPTION	
015-61622	Premium Ticonimum Wire 18-gauge	\$22.84

#### **Dental Alloy by Mesa distributed by ATI**

Cobalt based A alloy for prostheses, type 5

GROUP/ITEM #	DESCRIPTION	
110-40090025A	Magnum VIP-A Cobalt Alloy, 1 kg	\$310.00

### THE HISTORY OF POZZI ORTHODONTIC SILVER SOLDER Malibu



The Pozzi family has been manufacturing silver solders since 1934. From the melting of different alloys, including 65% of the purest silver, up to the making of rough bulk ingots, the quality control has been consistent, and the resulting product is remarkable. Of course, our solder does not have any cadmium or any other impurities. The rolling of the rough ingots to obtain square and round wires takes at least fifteen steps. The lamination is done with the most precise

gauged rolling mill in the industry. The high content of silver in our orthodontic solder means that very little material is needed in the soldering process, and the possibility of oxidation is nearly eliminated. Solders with a low silver content tend to oxidize heavily due to impurities. The combination of the ingredients in our silver solder and the high quality flux make this product the very best available in the world market.

### Silver Solder

#### Pozzi<sup>®</sup> Malibu<sup>™</sup> Orthodontic Laboratory Silver and Brazing Solder CADMIUM & NICKEL FREE

A cadmium-free alloy with free flowing and ductile characteristics. Easy to use. Manufactured in the U.S.A. for the Pozzi brand.

GROUP/ITEM #	DESCRIPTION		SPECIAL PRICE
031-32.50.001	1oz Spool	\$135.25	\$94.30
031-32.50.005	5oz Spool	\$676.18	\$471.15
031-32.50.050	50oz Spool	\$6,762.04	\$4,715.00

#### **Silver Solder A50N**

With the addition of nickel the A50N offers good strength and flow characteristics and works well on stainless steel and carbide.

GROUP/ITEM #	DESCRIPTION		SPECIAL PRICE
A50N	Pozzi Ortho Silver Solder .032 – 1oz	\$43.60	\$39.60
A50N-5BLK	Pozzi Ortho Silver Solder .32 – 5oz	\$201.83	\$196.20





Diameter	Melting Point	Silver Content
1/32" .79 mm)	1145 <sup>0</sup> F 618 <sup>0</sup> C	<b>56</b> %

Good strength and ductility on nickel alloys to help avoid stress corrosion cracking. Good color match on stainless steel. Excellent fluidity that promotes penetration for good strength and elongation

Offer Expires December 31, 2024–not valid with other offers www.americantooth.com | buy online: store.americantooth.com



## Pages 62-70

### **Flat Nose Pliers**



#### Wide Flat Nose Pliers

General utility pliers with squared beak. Heavy duty.

GROUP/ITEM # DESCRIPTION 007-8-601-001 Wide Flat Nose Pliers 51/2" (14 cm), serrated \$92.79 007-8-601-002 Wide Flat Nose Pliers 51/2" (14 cm), smooth \$93.25



### **Round Nose Pliers**

Same pliers as 007-8-601-022 with larger beaks.

GROUP/ITEM # DESCRIPTION 007-8-601-023 Round Nose Pliers 51/2" (14 cm), large tip \$92.80



GROUP/ITEM # DESCRIPTION 007-8-601-038

and Waldron.

**Triplex** 

Triplex 51/2" (14 cm) \$163.24



#### **Buy 3 Pliers or Instruments Same Type Get 1 FREE SAME TYPE**

#### Small Flat Nose Pliers

General utility pliers for detail work.

GROUP/ITEM # DESCRIPTION 007-8-601-006 Small Flat Nose Pliers 43/4" (12 cm), serrated \$88.08 Small Flat Nose Pliers 43/4" 007-8-601-007 (12 cm), smooth \$83.23

### **Round Nose Pliers**

Utility pliers specifically designed with flat-edged tips and round beaks on both sides.

GROUP/ITEM # DESCRIPTION 007-8-601-022 Round Nose Pliers 51/2' (14 cm), small tip \$90.89

#### Peeso

Band and crown stretching plus contouring pliers.

GROUP/ITEM # DESCRIPTION 007-8-601-051 Peeso 6" (15 cm) \$121.59

### **207272** Special Guarantee and Maintenence Guide

Pozzi Stainless Steel Pliers are the best quality instruments your money can buy! Made from the highest grade stainless steel which affords you maximum cutting and gripping capacity. These quality pliers being highly resistant to corrosion, are designed to stand up to frequent sterilization procedures, both hot and cold. Our pliers and instruments are guaranteed to be free of defects in workmanship and materials or it will be replaced at our descretion free of charge. This guarantee is for a period of two years unless stated otherwise.

Instrument Care: To insure maximum life expectancy for your instruments, please follow these recommended care instructions. When cold sterilizing solutions are used, never leave instruments soaking overnight; 30 minutes is sufficient. Remove instruments and towel dry them immediately. As a precaution, it is recommended to air dry the screw joints carefully with compressed air.


Caution: To assure the adjustable screw joints or Internal Pin Joints do not get rusty, DO NOT ALLOW THE INSTRUMENTS TO SELF DRY!

When autoclaving instruments the temperature should never exceed 200° C or 392°F. Once sterilizing is complete, always towel and air dry as recommended in the cold sterilizing procedure.

NOTE: We highly recommend applying a thin coat of silicon lubricant both to joints and tips at least once a week. WARNING: At time of shipping from our warehouse, this product is NOT sterilized. If this product is to be used intraorally it MUST be sterilized by an approved method prior to use.

Offer Expires December 31, 2024-not valid with other offers www.americantooth.com | buy online: store.americantooth.com



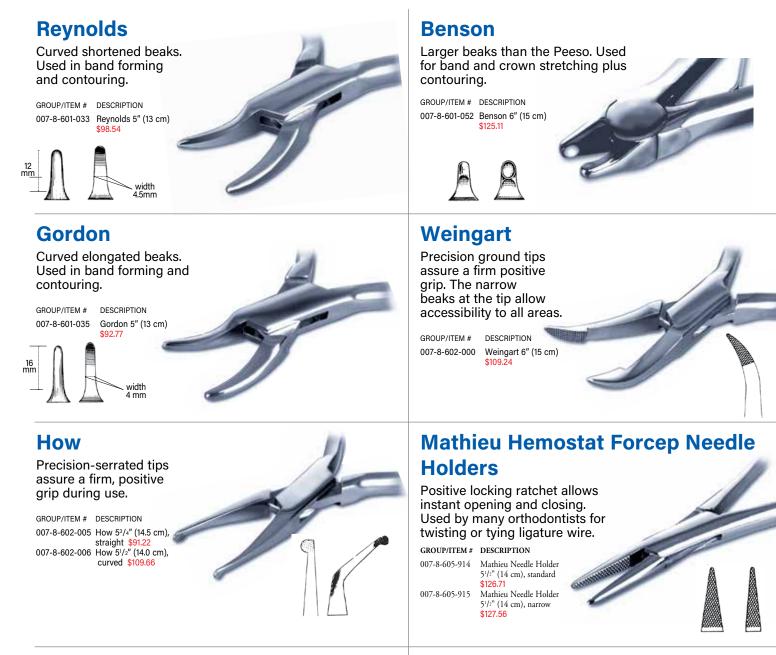

### **Broken Broaches Removing Pliers**

Elongated and tapered tips for removing elastics from inaccessible areas.

GROUP/ITEM # DESCRIPTION

007-8-601-101 Broken B Removir 5<sup>1</sup>/2" (14

DESCRIPTION
 Broken Broaches
 Removing Pliers
 5<sup>1</sup>/2" (14 cm), curved
 \$117.65



### **Anterior Band Removing Pliers**

Designed for the removal of anterior bands. Prevents danger of stretching or deforming the band. Comfortable for the patient without marring teeth.

GROUP/ITEM # DESCRIPTION 007-8-602-535 Anterior Band Removing Pliers 5" (13 cm) \$90.59



#### **Pozzi® Utility Pliers**

### **Posterior Band Removing Pliers**

**Nance/Adams Clasp Pliers** 

Extended jaw design allows the removal of posterior bands from inaccessible areas. Plastic tips are optional.

GROUP/ITEM # DESCRIPTION 007-8-602-536 Posterior Band Removing Pliers 5<sup>1</sup>/<sub>2</sub>" (13.5 cm) \$122.16

Allows the technician to

GROUP/ITEM #

007-8-601-034

GROUP/ITEM #

007-8-605-005

produce a nance or adams

clasp quickly and efficiently.

DESCRIPTION

Hemostat / Forcep

prevents crushing of the ligature.

DESCRIPTION

Mosquito Tip

Straight \$148.20

Hemostat / Forcep

**Mosquito Tip Straight** 

Modules sit securely in the tip. Hemostat

\$151.33

Nance/Adams Clasp Pliers



#### Johnson

Crimping and contouring pliers for bands, crowns and general use.

GROUP/ITEM # DESCRIPTION 007-8-601-028 Johnson 51/2" (14 cm) \$95.78



### **Anterior Direct Bond Bracket Remover**

Wedging the tips under the bracket brace, allows for easy and efficient removal. Firm grip assures ease in horizontal use.

GROUP/ITEM # DESCRIPTION 007-8-601-074 Anterior Direct Bond Bracket Remover \$114.16

# Hemostat / Forcep

Provides easy access to posterior

GROUP/ITEM # 007-8-630-8014

DESCRIPTION Hemostat / Forcep Mosquito Tip Curved \$137.21



# **Mosquito Tip Curved**

regions.

#### Hemostat / Forcep **Mosquito Hook Tip**

Hemostat has a recessed holding hook that prevents tearing of either the ligature or elastic.

GROUP/ITEM # 007-8-605-004

DESCRIPTION Hemostat / Forcep Mosquito Hook Tip \$154.29



### **Rose Torguing Pliers**


Rose Torauina Pliers wide with key .018 x .022. Designed to place torque in a specific section of archwire (single tooth) without distorting the remainder of the wire. The key is designed to torgue the section between holding beaks - the torqued section will not extend past the contact points of the pliers beaks and will engage easily without transfer of torgue to adjacent teeth. Labial torque and lingual torque can be placed next to each other with ease. Works with all arch wires.

GROUP/ITEM # 007-OM-805W DESCRIPTION Rose Torquing Pliers Wide with key \$126,79



Offer Expires December 31, 2024-not valid with other offers www.americantooth.com | buy online: store.americantooth.com CALL TODAY TO PLACE YOUR ORDER 805-487-9868

#### **Pozzi® Wire Benders**

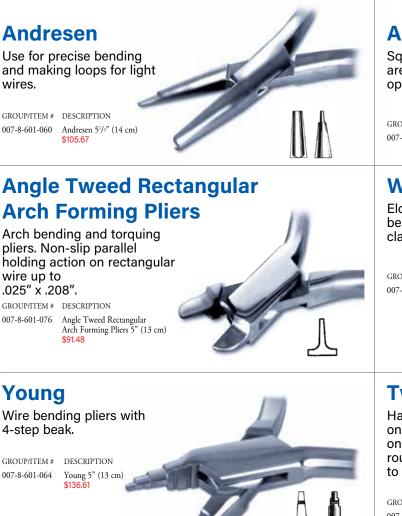

#### **Aderer Three Pronged Pliers**

Use for opening, closing and adjusting clasps.

GROUP/ITEM # DESCRIPTION 007-8-601-040 Aderer Three Pronged 4<sup>3</sup>/4" (12 cm) \$108.55



#### Adams

Squared tips. The jaws are parallel at 1 mm opening.

GROUP/ITEM # DESCRIPTION 007-8-601-043 Adams 5" (13 cm) \$93.36

#### **Wire & Clasp Bending Pliers**

Elongated beaks for precise bending of light wires and clasps.

GROUP/ITEM # DESCRIPTION 007-8-601-079 Wire and Clasp Bending Pliers 5<sup>1</sup>/<sup>2°</sup> (14 cm) **\$102.40** 



#### **Tweed Loop**

Has narrow concave beak on one side with three step on the other side for use on round or edgewise wire up to .020".

GROUP/ITEM # DESCRIPTION 007-8-601-070 Tweed Loop 5<sup>1</sup>/2" (14 cm) \$106.41

### Aderer

Same pliers as 007-8-601-040 with smaller and shorter tips.

GROUP/ITEM # DESCRIPTION 007-8-601-048 Aderer Three Pronged 4<sup>1</sup>/-<sup>2</sup> (11.5 cm) \$94.99



#### Pozzi 3-Jaw Plier

For bending of Rapid Expander Arms, for wires up to 1.5 mm.

GROUP/ITEM # DESCRIPTION 019-AP1034-00 3-Jaw Round Plier (Sand Finish) \$51.75





### **Pozzi® Wire Benders**

#### Aderer

Arch bending pliers. Two types - 2 groove and 3 groove.

GROUP/ITEM # DESCRIPTION 007-8-601-012 Aderer 5<sup>1</sup>/4" (14.5 cm) 3 grooves \$98,49 TA

007-8-601-024 Aderer 5<sup>1</sup>/2" (14 cm) 2 grooves \$99.02

### Aderer Lingual Bar Heavy Face Bow

Three-pronged box joint pliers designed for the heavy bending of lingual bars, face bows and clasps. Heavy duty.

GROUP/ITEM # DES 007-8-601-041 Ade (He

DESCRIPTION Aderer Lingual Bar (Heavy Face Bow) \$101.33



### **Bird Beak Pliers**

Both round and flat beaks permit the bending of rectangular and round wire up to .030"

GROUP/ITEM # DESCRIPTION 007-8-601-077 Bird Beak 5" (13 cm) \$111.73



### Lite Wire

Pliers designed for forming and bending up to .020" wire. Also for making and closing loops.

GROUP/ITEM # DESCRIPTION 007-8-601-078 Lite Wire 5" (13 cm) \$105.05



#### Schwarz

Light wire, small optical pliers with concave beaks.

GROUP/ITEM # DESCRIPTION 007-8-601-067 Schwarz 5" (13 cm) \$93.51



### Waldron

Arch and spring bending pliers.

GROUP/ITEM # DESCRIPTION 007-8-601-030 Waldron 5" (13 cm) \$93.43



### Jarabak

For precise bending of light wires. The three grooves ensure accurate bending and closing of loops.

GROUP/ITEM # DESCRIPTION 007-8-601-063 Jarabak 5<sup>1</sup>/2" (14 cm) \$101.95



### **Hollow Chop**

For contouring round arch wires. Smooth beaks. GROUP/ITEM # DESCRIPTION 007-8-601-069 Hollow Chop 5<sup>1</sup>/\*" (14 cm) \$97.39



### **Optical Pliers**

Precision concave beaks for clasp and wire bending. Can be used on both round and edgewise wire to bend to exacting measurements.

GROUP/ITEM # DESCRIPTION 007-8-601-071 Optical 4¾" (13 cm) \$91.36



### **Wire Bending Pliers**

Precision round and concave beaks for clasp and wire bending. Can be used on both round and edgewise wire up to .040.

GROUP/ITEM # DESCRIPTION 007-440 Wire Bending Pliers \$134.41

> CALL TODAY TO PLACE YOUR ORDER 805-487-9868

**Pozzi® Wire Benders** 

### **Clasp Bending Pliers**

007-442

GROUP/ITEM # DESCRIPTION **Clasp Bending Pliers** \$136.34



### **Heavy Duty Wire & Plate Scissors**

GROUP/ITEM # DESCRIPTION Heavy Duty Wire and Plate Scissors 017-437

\$83.19



### Lab Cutting Pliers

Designed to cut lab wires up to 0.040".

GROUP/ITEM # DESCRIPTION 007-8-603-091H Lab Cutting Pliers \$309.36



### **Wire Cutting Pliers**

Designed to cut lab wire up to 0.040".

GROUP/ITEM # DESCRIPTION 007-8-603-092 Wire Cutting Pliers \$110.23



### **Ligature Wire Cutter**

Designed to clip metal ligatures.

GROUP/ITEM # DESCRIPTION 007-8-603-093 Ligature Wire Cutter \$107.50



### **Heavy Duty Waldsachs**

General forming and bending pliers. Heavy duty.

GROUP/ITEM NO. DESCRIPTION 007-443 Cutting Pliers, Waldsachs \$139.21



### **Scissors And Cutters**

Universal collar and crown scissors. For cutting bands and soft wire.

GROUP/ITEM # DESCRIPTION 007-8-604-046 Crown Scissors, Smooth Edge 43/4" (12 cm), Straight \$89.52

007-8-604-047 Crown Scissors, Smooth Edge 43/4" (12 cm), Curved \$97.39



### **Pin & Ligature** Wire Cutter

Special carbide steel blades assure accurate, dependable cutting edges. Perfect for cutting lock and key pins. Stainless steel. Cuts soft wire up to .010". Do not use on hard wire.

GROUP/ITEM # DESCRIPTION 013-19790 Pin and Ligature Wire Cutter \$137.95



#### **Hemostat/Forcep Mosquito Hook tip**

Hemostat has a recessed holding hook that prevents tearing of either the ligature or elastic.

GROUP/ITEM # DESCRIPTION 007-8-605-004 Hemostat / Forcep Mosquito Hook Ťip \$154.28



### **Heavy Duty Lab Pliers**

Designed to cut lab wire up to 0.032".

GROUP/ITEM # DESCRIPTION 007-8-603-096 Heavy Duty Lab Pliers \$135.94



CALL TODAY TO PLACE YOUR ORDER 805-487-9868

Offer Expires December 31, 2024-not valid with other offers www.americantooth.com | buy online: store.americantooth.com

#### Pozzi<sup>®</sup> Classika

We wanted to manufacture a pliers like the ones you used to be able to buy in the old days - without having to mortgage the farm!

We feel that once you use one of these pliers you will know why it is a real classic!

All of the "Classika" instruments are made of corrosion resistant stainless steel. All "Classika" cutting instruments have tool steel inserts. All "Classika" cutting instruments have diamond-honed cutting edges and have a hard chrome finish. They withstand all accepted

methods of sterilization.

All "Classika" instruments have a balanced assembly wobble-free sliding lock joint. No more shim joints to deal with. All come with extended "sure grip" handles and are guaranteed to be free from manufacturing defects. Any instrument that proves defective in workmanship and materials will be repaired or replaced at our discretion free of charge. All instruments carry a two year warranty.

### **Black Hard Wire Cutter**

Specially designed for cutting hard wires. The cutting angle is 15 degrees and the cutting edges have been diamondhoned to allow hard wire cutting. Maximum cutting capacity is .20 diameter in hard wire.

GROUP/ITEM # DESCRIPTION 040-40-0015 Black hard wire cutter \$181.19



# Flush Cut Distal End Cutter with Safety Hold

The particular design of this cutter allows distal end arch wires to be cut in the patient's mouth safely. It cuts wire and then "safety holds" the loose wires.

Maximum cutting capacity: rectangular wire up to .022x.028"(0,56x0,71mm); round wire up to .020" (0,51mm); twist wire up to .016" (0,41mm) in diameter.

GROUP/ITEM # DESCRIPTION 040-40-0017R Distal end cutter \$181.19

### Miniature Pin & Ligature Cutter

This standard pin cutter is non-angled and designed to cut soft wires. This pliers is very popular for chairside work. Maximum cutting capacity is .015 in soft wire.

GROUP/ITEM # DESCRIPTION 040-40-0016 Miniature pin and ligature cutter \$181.19



#### Micro-Miniature Light Wire Cutter with Stop Screw

This pliers is designed as a miniature ligature cutter that will cut soft wire. It easily reaches those hard-to-reach areas in the mouth. It has a factory adjusted stop screw that prevents tip damage due to overscrewing the cutter during routine procedures. Maximum cutting capacity is soft wire up to .012 in diameter.

GROUP/ITEM # DESCRIPTION 040-40-0018 Micro-Miniature cutter \$18119 15° Pin & Fine Wire Cutter

The tip is angled 15 degrees for better field vision. The tooled steel inserted tips have backing ledges for greater strength. Maximum cutting capacity is .015 diameter for both soft and hard wire.

CALL TODAY TO PLACE YOUR ORDER

805-487-9868

GROUP/ITEM # DESCRIPTION 040-40-0014 15 degree pin and fine wire cutter \$181.19

## Offer Expires December 31, 2024-not valid with other offers

#### **Stainless Steel Crown Cutting Tools**

#### **Small Straight Scissors**

Use for trimming anterior crowns.

GROUP/ITEM # 040-40-0028 \$140.28



**Tube Crimping Pliers** 

Used to crimp the tube around a space maintainer wire. It can also be used for placing stops and tubing on arch wires.

GROUP/ITEM # 040-40-0033 DESCRIPTION Tube Crimping Pliers \$246.93



#### **Stainless Steel Crown Pliers**

**Small Crown Crimping Pliers** This tool is handy for anterior crowns and smaller sizes in the primary molars.

GROUP/ITEM # 040-40-0031 \$246.93



**Regular Crown Crimping Pliers** This is the standard in the industry for all posterior crowns.

GROUP/ITEM # 040-40-0030 \$246.93 **Crown Contouring Pliers** For adding contour to the crown for a better fit.

GROUP/ITEM # 040-40-0032 \$246.93



CALL TODAY TO PLACE YOUR ORDER

805-487-9868

#### **Mini-Distal End Cutter**

Smaller than standard. Cuts and holds in one operation.

**GROUP/ITEM #** 007-8-603-090H

DESCRIPTION Mini-Distal End Cutter



### Posterior Direct Bond Bracket Remover

Direction provides easy access to confined posterior areas. Easy vertical removal.



### **Boynton-style Elastomeric Pliers**

Features a smaller head with narrower arms for easier access to the posterior. An ideal alternative multiple use tool for the Crown & Bridge and Ceramics department.



### **V-Stop Bend Pliers**

Forms a 1 mm high V bend automatically. Used to shorten arch wires or to place a strategic stop in the arch wire, such as the midline, which avoids any unwarranted displacement of the archwire. Used by laboratories for wrought wire modifications in transitional partials.

For wires up to .025".

GROUP/ITEM # 007-8-601-045 \$121.36

### **Side Cutter Pliers - Heavy Duty**

Heavy duty pliers designed to cut wire up to .040"

GROUP/ITEM # 007-8-603-096H \$221.96



#### **Side Cutter Pliers**

Will cut up to .020" round wire. Features dual reciprocal springs for a lighter touch and less hand fatigue.

GROUP/ITEM # 007-8-603-098H \$112.55



### Mathieu Needle Holder

Positive locking ratchet allows instant opening and closing. Used by many orthodontists for twisting or tying ligature wire. An excellent general use instrument for Crown & Bridge and Ceramic technicians worldwide.



### **Prolong Flat Nose Pliers**

#### Serrated

Similar to a standard Mathieu except the tips are ultra fine to facilitate application of elastomeric ligature ties.

GROUP/ITEM # 007-8-605-125 \$129.52



### Mathieu Needle Holder

Short Beak

General detail serrated utility pliers for detail work 6" - 15cm

GROUP/ITEM # 007-8-601-015 \$88.60



### Pozzi Special Guarantee

Pozzi pliers and instruments are guaranteed for a period of two years from date of purchase to be free from defects in workmanship and/or materials. Any product that proves defective in workmanship and materials will be repaired or replaced, at our discretion, free of charge.



### **TCS Dental Supply Products Distributed by ATI**

#### SPECIALIZED DENTURE LABORATORY SUPPLIES

## Unbreakable

#### Thermoplastic Comfort Systems (**tcs**<sup>®</sup>)

tcs unbreakable partials are fabricated from a strong and durable biocompatible nylon thermoplastic. tcs unbreakable may be used whenever you are considering a removable partial denture for your patient. tcs unbreakable gives you the option of offering your patients a restoration with flexiblity, strength, and aesthetics.

#### Some benefits included:

- Color stable
- Ideal for combo partials
- Lifetime warranty against breakage Will not warp or become brittle
- Made in the U.S.A.

# Karadent

Karadent<sup>™</sup> was introduced to the market in 2016.

It is a semi-rigid microcrystalline polymer for full dentures, semi-flexible partial dentures, clear clasps, frameworks & nightguards.

- No volume shrinkage exceptional fit
- Minimal water absorption
- An upgrade from acrylic dentures
- Quick and easy to finish and polish (similar to acrylic)
- · Easy to repair, rebase, and reline with acrylic (laboratory procedure using fusing liquid)
- Monomer free and hypoallergenic
- BPA free
- Available in 2 translucent shades

# iFlex

iFlex<sup>™</sup> is a polyolefin thermoplastic used for the fabrication of flexible partial dentures. Ideal for technicians new to flexible partials Clinically unbreakable Easy to finish and polish Stain-free High translucency Monomer free and hypoallergenic BPA free. Available in 5 translucent shades Cartridges: iFlex™

		Made in	
SMALL - 5 per Pkg	\$32.55	004	
Group/Item#	Description	SPECIAL PRICE	10766
099-4002-01/S Shade*	<b>Removable Partial</b>		
	material 25 mm	\$28.00	
MEDIUM - 5 per Pkg	\$43.93		
Group/Item#		SPECIAL PRICE	WH I
099-4000-01/M Shade*	<b>Removable Partial</b>		1111
	material 25 mm	\$38.00	ΠH
LARGE - 5 per Pkg	53.25		
Group/Item#		SPECIAL PRICE	HI
099-4005-01/L Shade*	Removable Partial material 25 mm	\$46.00	
*Shadoe: Light Pink-1 St	tandard Pink-2 Light	/Dark	111-1

\*Shades: Light Pink=1, Standard Pink=2, Light/Dark Pink=3, Dark Pink=4, Natural=N



LARGE - 5 per Pkg...... \$53.25

Group/Item#

5005-01/L Shade\* **Removable Partial** material 25 mm

\*Shades: Light Pink=1, Standard Pink=2, Light/ Dark Pink=3, Dark Pink=4, Natural=N



Offer Expires December 31, 2024-not valid with other offers www.americantooth.com | buy online: store.americantooth.com CALL TODAY TO PLACE YOUR ORDER 805-487-9868

### **TCS Dental Supply Products Distributed by ATI**

#### SPECIALIZED DENTURE LABORATORY SUPPLIES



## **Pour Flask**

The Pour Flasks are used for the acrylic pour technique and to duplicate models. They are manufactured in light-weight aluminum. Easy to assemble, disassemble and clean.

Group/Item# Description 099-3650-01 Pour Flask



## **Econo Cutter Mandrel**

Special mandrel required for the 4" diameter Econo Cutter



Group/Item# 080-4204-01

**Description** Econo Cutter Mandrel

\$11.20

## **Econo Cutter X-Lrg.**

\$46.92

The 4" Diameter Econo Cutter is used to cut off the sprues from the tcs<sup>®</sup> restorations.



\$11.20

Group/Item# 080-4203-01 **Description** Econo Cutter X-Lrg.

## Knock Down Wheel Chuck

Mandrel used with the White Knock Down Wheel and the Rubberized Abrasives (wheel 304F and 306F).

Group/Item# 099-4212-01 Description Knock-Down Wheel Chuck

\$23.36

## White Knock-Down Wheel

Ideal for tcs<sup>®</sup> restoration as one of the first steps in finishing. Use to cut the excess bulk and trim partial to the desired length.

Group/Item# 080-4211-01 Description
White Knock-Down Wheel 3 x 3/8 \$16.95

## **Knife Blades**

Replacement Blades for the tcs<sup>®</sup> Carving Knife.

Group/Item# 099-4202-01

# Description Knife Blades - 12 Pack

Pack \$25.42

## **Brown Tripoli**

Designed for high shine of tcs° unbreakable.

Group/Item# 099-4206-01

Description Brown Tripoli - 2 lb \$26.67



An advanced formula for high-speed polishing of acrylics and thermoplastics

**Group/Item#** 099-4298-01 099-4298-25 Description Polycril 1 lb \$4.38 Polycril 25 lb \$88.48



### **TCS Dental Supply Products Distributed by ATI**

SPECIALIZED DENTURE LABORATORY SUPPLIES

## **Carving Knife**

This knife is excellent for removing extended fibers that sometimes appear during processing of tcs<sup>®</sup> partials. It is also used to define the necklines of the tcs® restorations when necessary.



\$32.58

**GROUP/ITEM #** 099-4201-01

DESCRIPTION tcs® Carving Knife Rag Buff

High-quality and long-lasting rag buff stitched only in the center for the extra softness needed for polishing tcs® restorations. The rag buff is used dry with pumice or pumice substitute for the high-shine process. Diameter 4" x 40 ply



#### GROUP/ITEM # DESCRIPTION

099-4207-01 4"x40 ply Leather Ctr Buff - Hi Shine \$5.17

## **Rubberized Abrasives**

Long-lasting wheel for smoothing tcs® restorations during the finishing stage.



GROUP/ITEM #	DESCRIPTION	
099-4304-01	Rubber Wheel 304-F	\$20.83
099-4306-01	Rubber Wheel 306-F	\$24.17

## **Rubber Wheels**

The Small Brown Rubber Wheels are used as the last step prior to polishing tcs<sup>®</sup> restorations.




**GROUP/ITEM #** DESCRIPTION Rubber Wheel #723 099-4503-01 \$27.92

## **Green High Shine**

Designed for high shine of tcs<sup>®</sup> unbreakable.



**GROUP/ITEM #** DESCRIPTION \$23.80 099-4209-02 Green High Shine - 1lb

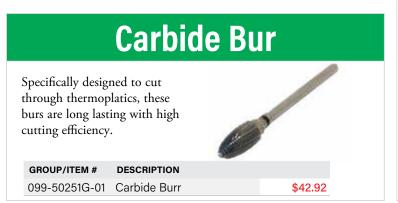

Offer Expires December 31, 2024-not valid with other offers www.americantooth.com | buy online: store.americantooth.com

## **Pink Grinding Wheels**

The Pink Grinding Wheel is an excellent and economical product for the finishing stage of tcs° restorations.



GROUP/ITEM #	DESCRIPTION	
099-4205-01	Pink Grinding Wheels	\$28.00



#### BUY 3 Wax & Modelling Carvers **ASA DENIAL** MADE IN ITALY **1 FREE GET** WAX AND MODELLING CARVERS Group / Item # 017-5100-1 Graded Handle Hylin Carver 0.0.0.0.18.18.18.18.18.18.18.1 \$13.92 ASA DENTAL 5100-1 C€ 017-5100 Hylin \$13.92 ASA DENIAL 5100 CE 017-5101 -0 Hylin-Le Cron 5101 C€ \$13.92 017-5102 Hylin-Ceramica \$13.92 ASA DENTAL 5102 €€ \*\*\*\*\*\*\*\*\*\*\*\*\*\*\* 017-5103 Hylin-ASA \$13.92 ASA DENTAL 5103 C€ -----017-5104 Wax-Model Carver \$13.92 ASA DENTAL 5104 CE 017-5105 Le Cron \$13.92 ASA DENTAL 5105 CC ANNANANANANAN 017-5106 Le Cron \$13.92 COLUMN COOL ASA DENTAL 5106 CC 017-5107 αŝ Le Cron \$13.92 ASA DENIAL 5107 CC 017-5108 Le Cron \$13.92 ASA DENTAL 5108 CE. ..... 017-5109 Zahle \$13.92 017-5110 Zahle ------\$13.92 017-5111 Zahle \$13.92 017-5113 Gritman ASA DENIAL 5113 C€ \$13.92 017-5114 Beale \$13.92 ASA DENTAL 5114 CE THE REPORT OF THE PARTY OF THE 017-5115 Roach RSR DENTRL 5115 CC \$13.92 017-5116 Vehe Reporter of \$13.92 ASA DENTAL 5116 CC

74

### ASA Dental Products Distributed by ATI



### SUITABLE "GALLONI" CRUCIBLES FOR "CMP INDUSTRIES LLC" CASTING MACHINES



#### CERAMIC CRUCIBLE

For Chrome-Cobalt-Molybdenum Alloy. Suitable for last models: Modular X and Modular 20 (not in use on other models).

GROUP/ITEM #DESCRIPTION108-2010MP4Ceramic Crucible\$34.00



#### CERAMIC SILICON NITRIDE CRUCIBLE WITH HIGH RESISTANCE

For Chrome-Cobalt-Molybdenum and Palladium Alloys. Suitable for last models: Modular X and Modular 20 (not in use on other models).

GROUP/ITEM #DESCRIPTION108-2010MP5Ceramic Silicon Nitride Crucible with High Resistance\$30.97



#### CERAMIC CRUCIBLE

For Chrome-Cobalt-Molybdenum Alloy. Suitable for all previous models. Suitable for all old machine models.

GROUP/ITEM #DESCRIPTION108-4114Ceramic Crucible\$15.20



#### CERAMIC SILICON NITRIDE CRUCIBLE WITH HIGH RESISTANCE

For Chrome-Cobalt-Molybdenum and Palladium Alloys. Suitable for all old machine models.

**GROUP/ITEM #DESCRIPTION**108-4115Ceramic Silicon Nitride Crucible with High Resistance\$22.90



GRAPHITE CRUCIBLE COMPLETE WITH CEMENTED CERAMIC CONTAINER AND LID COVER

For Precious Alloys. Suitable for all old machine models.

GROUP/ITEM #	DESCRIPTION	
108-4116	Crucible complete with cemented ceramic container and lid cover	\$58.83





For Chrome-Cobalt-Molybdenum Alloy. Suitable for all old machine models. Smaller Capacity (45 grams).

GROUP/ITEM #DESCRIPTION108-3606Ceramic Crucible\$10.38

We also have parts for CMP Industries LLC casting machines please inquire with ATI.



# **Shining Box** DENTURE POLISHING STATION



SHINING BOX is a complete cabinet developed for grinding and polishing operations under full safety conditions on metals and other materials currently used in the dental field.

The hood has been designed and constructed with the aim to protect the user against any projection of fluids and dust particles and to extract and filter the residues. Configuration of the polishing unit provides for dual speed settings: at low medium speed (1.400 rpm) when removing particles on metals without overheating in excess the surfaces and medium-high speed (2.800 rpm) for excellent polishing effect. The conic tips can accommodate abrasive wheels, brushes and fabric discs.

The highest safety during operation is assured by the adoption of an electromagnetic switch which prevents the unit from starting up accidentally after a voltage drop or lack of current. In addition, we recommend to wear specific individual safety protection devices as goggles, gloves and mask to preserve your eyes from any slag rebounds, your hands from touching accidentally the parts in motion and your respiratory tract from inhaling dust particles.

#### **TECHNICAL FEATURES**

Speed: 1,400/ 2,800 rpm Absorption: 300W Suction power: 500 mc/h Noise level: >70 dBA Lighting: LED Overall dimensions: mm 730x 500x 580 H Total Weight: 49 kg

GROUP/ITEM #	DESCRIPTION	RETAIL PRICE	SPECIAL PRICE
109-A5250	Denture Polishing Station 110V	\$6,508.68	\$5,500.00



#### **Dentalfarm products distributed by ATI**

# BASE 3



The **BASE 3** is the most popular and successful model of the entire and wide range of sandblasters that Dentalfarm has ever manufactured. It is a complete blasting center consisting of:

- one working chamber for the coarse aggressive corundum used to remove investment, propelled by a powerful, recycling projector.
- two disposable microblasting EASY tanks for use with aluminum oxide or polishing microbeads.

The standard range already incorporates the necessary instrumentation and controls, recently implemented by a new LED white lighting system.

#### **TECHNICAL FEATURES**

**Recycling blasting system:** devesting projector fitted with 1x0 3.0 mm nozzle.

Microblasting tanks: two (EASY version)

**Available nozzles:** 0.8 mm for grainsize from 75 to 95µ; 1.2 mm from 95 to 110µ; 0.5 mm for 50µ; 2.0 mm for 250µ

Optional connection to an electric extractor (PRD-3 or PRD-3 SHAKE) or to the Dentalfarm WAFI.S. solution Working pressure: from 2.5 -6.0 BAR (devesting) from 1.5 -5.0 BAR (micro blasting)

Overall dimensions: 415x400x415mm

Blasting chamber volume: 22 I

Weight: 11.8 kg

GROUP/ITEM	# DESCRIPTION	RETAIL PRICE	SPECIAL PRICE
109-A10723	All-in-one blasting center	\$2,286.68	\$1,652.13

# MICRA 2



Microblaster standard fitted with two microblasting tanks for aluminum oxide or micro glassbeads for final polishing. The working chamber provides for a collection compartment for disposable abrasive products, separated by a filtering grate. The front large viewing section is fully transparent allowing for greater visibility inside the blasting chamber. Totally integrated unit which accomodates the necessary instrumentation and controls providing for ideal operation, now implemented by a new LED lighting system.

MICRA 2 offers the feasibility of two different connections for the suction systems: an upper intake to collect dust particles only, and a lower intake to collect both dust and used abrasive material. The MICRA EVO and EVO D.O.S. versions can be upgraded to up to 4 EASY and D.O.S. tanks.

#### **TECHNICAL FEATURES**

Microblasting tanks: 2 (EASY version)

**Available nozzles:** 0.8 mm for grainsize from 75 to 95µ; 1.2 mm from 95 to 110µ; 0.5 mm for 50µ; 2.0 mm for 250µ

**Operation:** electric foot-control and selector switch

Hand sleeves: high sensitivity latex gloves

Optional connection to an electric extractor (PRO-3 or PRO-3 SHAKE) or to the Dentalfarm W.A.FI.S. solution

- Working pressure: 1.5 5.0 BAR
- Overall dimensions: 410x385x415mm

Blasting chamber volume: 24 |

Weight: 11 kg

GROUP/ITEM #	DESCRIPTION	RETAIL PRICE	SPECIAL PRICE
109-A10732	Finishing microblaster	\$2,180.22	\$1,575.21

## **MINIEXPORT**



Vibrator for plaster and investment, featuring a rectangular plate. The vibration is perfectly vertical without any dispersion thanks to 4 shock absorbers and special large suction feet providing uniform vibration. Intensity and width of movement of the plate are controlled by a circuit board fitted with a potentiometer.

All components are soaked in resin in order to assure the highest stability and insulation. The rubber plate is easily removable for cleaning. MINIEXPORT is ideal for small dental labs.

Plate size: 240x150x105mm

It accommodates 2 large cylinders (9x).

#### Weight: 5,5 kg

GROUP/ITEM #	DESCRIPTION	RETAIL PRICE	SPECIAL PRICE
109-A0125	Electronic regulation capacity 2 chrome rings	\$ 628.23	\$453.90

## **EXPORT 2**



Vibrator for plaster and investment, featuring a rectangular plate. Four corner post construction provides uniform vibration to all points of the tray allowing full use of the table area. Perfect stability is assured by special suction feet. Intensity and width of movement of the plate are controlled by a circuit board operated by a potentiometer. All components are soaked in resin in order to assure the highest stability and perfect insulation. The rubber plate is easily removable for cleaning. EXPORT 2 is ideal for large workload dental labs.

Plate size: 365x265x110 mm

It accommodates 6 large cylinders (9x).

#### Weight: 10,5 kg

GROUP/ITEM #	DESCRIPTION	RETAIL PRICE	SPECIAL PRICE
109-A0127	Electronic regulation capacity 6 chrome rings	\$ 865.17	\$625.08

## MAXIEXPORT



Vibrator for plaster and investment, featuring a very large rectangular plate. The vibration is perfectly vertical and without any dispersion thanks to 3 combined electromagnetic units and 8 shock absorbers assuring uniform vibration all over the plate. Large suction feet provide for stable operation. Intensity and width of movement of the plate are controlled by a circuit board operated by a potentiometer. All components are soaked in resin in order to assure the highest stability and perfect insulation. The rubber plate can be easily removed for cleaning purposes. MAXIEXPORT is ideal for large model casting dental labs.

Plate size: 600x400x130 mm

It accommodates 18 large cylinders (9x).

#### Weight: 24,7 kg

GROUP/ITEM #	DESCRIPTION	RETAIL PRICE	SPECIAL PRICE
109-A0129	Electronic regulation capacity 18 chrome rings	\$ 1,630.20	\$1,177.82

#### **PHOTO-FLASK**



#### PHOTOPOL PRO



#### IDEA



#### **IDEA PRO**



#### **Dentalfarm products distributed by ATI**

Transparent and made of removable sections, PHOTO-FLASK is a practical device designed to handle the light-curing composites without any modification of the standard operating working procedure. PHOTO-FLASK is made by transparent elements ensuring a better exposure to UV rays. No reinforcement metallic plates are needed.

GROUP/ITEM #	DESCRIPTION	RETAIL PRICE	SPECIAL PRICE
109-PHF	Flask for light-curing units	\$ 359.38	\$259.65

Photopol Pro is a a light-curing unite which best suits with most composites sensitive to UV light. 3 extremely powerful circuits and 8 LED UV at different frequency assure the overall radiation. Both the working chamber and the horizontal plate are coated with highly reflecting material assuring a mirror effect with no dispersion of ultraviolet rays. The support plate is steadily rotating and it has been conveniently assembled on a sliding drawer allowing to position and to pick up easily the elements to be treated.

#### **TECHNICAL FEATURES**

Two circuits of 8 LED each covering the higher UV spectrum (above 410nm-blue)

One circuit of 8 LED covering the lower UV in stead spectrum (under 410mm - violet) Rotatin

Adjustable irradiation power to avoid any shock to the material

Ideally suited to perform light-curing to any composite available in the dental field.

Useful "Step Mode" to dry out and fix different layers (operation lasting few seconds, repeated in steady succession)

Rotating turntable working under continuous cycle to assure the all-over exposure of the elements restricting the shade effect.

Plate size: 100x100 mm

Dimensions: 280x350x300 mm

Weight: 8.5 kg

GROUP/ITEM #	DESCRIPTION	RETAIL PRICE	SPECIAL PRICE
109-A5407	UV light-curing unit for composites	\$ 3,127.90	\$2,259.91

IDEA represents the utmost achievement in developing the design of the new range of sandblasters. This high-point combines the reliability of a brand like Dentalfarm, the ergonomics of the design and the modularity of the models available. The particular structure of the machine facilitates a comfortable position and assures the best angle of impact for the abrasive jet even when the user is in a seated position. The wide visibility over the entire workspace and the powerful lighting of natural color allows to perform very accurate finishing operations. The unique pointed jet will grant the maximum efficiency whilst limiting product consumption.

The possibility of adopting the heated D.O.S. (Dry Oxide System) tanks represents the best solution to problems originated by insufficient smoothness, especially for very fine abrasive powders. The sandblaster can be customized in both the number and pattern of containers (from two to four D.O.S or EASY tanks) and the type of dust extraction system (external electric or built-in pneumatic).

GROUP/ITEM #	DESCRIPTION	RETAIL PRICE	SPECIAL PRICE
109-A1077	Finishing microblaster	\$2,981.48	\$2,154.12

IDEA PRO is the response to the demand of every qualified dental laboratory. A multi-function solution that incorporates both devesting and finishing in the same unit featuring a design which avails of the latest technical solutions.

The entire working process embodied in a single unit: a preliminary stage of fast and powerful devesting and a second stage of careful and precise finishing.

The rationalization of space and the location of the tools allows to perform easily all blasting operations in the same working chamber.

This hybrid version has a very powerful jet that provides the recycled corumdum for devesting and a variable number of tanks for the finishing micro abrasives. The working chamber is equipped with LED lighting and separated from the collection tank by a fixed grid that prevents debris from mixing with the stored sand.

GROUP/ITEM #	DESCRIPTION	RETAIL PRICE	SPECIAL PRICE
109-A1076	All-in-one blasting center	\$3,111.92	\$2,248.36

#### DENTUREMART

A QUANTITY DISCOUNT PROMOTIONAL CATALOG VOLUME 2 • 2023-24 • ISSUE #29

#### Published by: American Tooth Industries 1200 Stellar Drive • Oxnard, CA 93033-2404 USA www.americantooth.com

#### **CALL TODAY TO PLACE YOUR ORDER** 805-487-9868

Buy Online: store.americantooth.com Email: info@americantooth.com Fax (805) 483-8482 Products are guaranteed. No other discounts apply.

PRSRT STD MAIL U.S. POSTAGE PAID OXNARD, CA PERMIT NO. 539

#### **AMERICAN TOOTH INDUSTRIES MISSION STATEMENT & QUALITY POLICY**

Our family of brands-Justi, Pozzi, Major, and Regal-carries an old motto of our company, "Simple inventions last forever." Using the tools of their forebears, the human eye and hand, our craftsmen carry on the American Tooth Industries' tradition of excellence. We are delivering quality products at a fair market value, and we shall continue to do so.

## BUY AMERICAN MADE



American Tooth Industries is now the ONLY Tooth Manufacturer in the USA

Justi has been proudly manufacturing Hardened Plastic Teeth for over 159 years in the US in an FDA Registered facility and continues to do so today. We have one of the largest tooth inventories in North America with over 6 million sets in stock.



PavPa